

**Complementary Technical Document**

# **UptiBlue Viable Cell Counting Kit**

## **Indications for use**

The UptiBlue Viable Cell Counting Assay is designed to measure quantitatively the proliferation of various human and animal cell lines, bacteria and fungi. The bioassay may also be used to establish relative cytotoxicity of agents within various chemical classes. The toxicologist can establish baseline data for predicting the toxicity of related novel agents by comparing such baseline data with known *in-vivo* toxicity.

The assay is simple to perform since the indicator is water soluble, thus eliminating the washing/fixing and extraction steps required in other commonly used cell proliferation assays.

## **Product Description**

The UptiBlue Viable Cell Counting Assay incorporates a fluorometric/colorimetric growth indicator based on detection of metabolic activity. Specifically, the system incorporates an oxidation-reduction (REDOX) indicator that both fluoresces and changes color in response to chemical reduction of growth medium resulting from cell growth.

The specific (fluorometric/colorimetric) REDOX indicator incorporated into UptiBlue has been carefully selected because of several properties. First, the REDOX indicator exhibits both fluorescence and colorimetric change in the appropriate oxidation reduction range relating to cellular metabolic reduction. Second, the REDOX indicator is demonstrated to be minimally toxic to living cells. Third, the REDOX indicator produces a clear, stable distinct change which is easy to interpret. The REDOX indicator has no current or past indication of carcinogenic capacity.

As cells being tested grow, innate metabolic activity results in a chemical reduction of UptiBlue . Continued growth maintains a reduced environment while inhibition of growth maintains an oxidized environment. Reduction related to growth causes the REDOX indicator to change from oxidized (non-fluorescent, blue) form to reduced (fluorescent, red) form.

Experiments suggest that reduction of UptiBlue requires uptake by the cells. To test this hypothesis, we grew A549 cells to confluency in T25 flasks using RPMI 1640. The media was then removed from 2 flasks containing cells and replaced with fresh media. Fresh media was also added to a sterile flask containing no cells to serve as a negative control. All flasks were then re-incubated at 37°C, 5% CO<sub>2</sub> for 4 hours.

At the end of the 4 hour incubation, UptiBlue was added to each flask. There was no immediate color change in any flasks upon addition. In one of the flasks containing cells, the media was left in contact with the cell layer, while the other T flask was turned over so that the media was not in contact with the cell layer. All flasks were then incubated for 1 hour at 37° and rechecked for color change.

If UptiBlue reduction occurred simply from the reduction of the external medium, we would expect the flask in which the media was in contact with the cells and the flask in which media was no longer in contact with the cells to exhibit the same amount of reduction. This was not the case. The flask where the media was not in contact with the monolayer following addition of UptiBlue displayed no color change from the blue of the negative control flask. The flask where the cells were in contact with the monolayer was very pink, indicating a higher percentage of reduction. This seems to indicate uptake by the cells is required for reduction of UptiBlue .

We have analyzed the possible interaction of UptiBlue in cellular respiration by looking at the relative redox potential of the various components of the biological respiration chain. These are presented in Table 1, which includes the oxidation-reduction potentials for UptiBlue and MTT and where they interact with the system.



**Table 1. Oxidation-reduction potentials in the electron-transport system.**

Half-reaction		E' <sub>o</sub> (mV) pH 7.0 25°C
O <sub>2</sub> + 4H <sup>+</sup> + 4e <sup>-</sup>	2H <sub>2</sub> O	+820
UptiBlue <sub>ox</sub> + 2H + 2e <sup>-</sup>	UptiBlue <sub>RED</sub>	+380
cytochromes <sub>ox</sub> + 1e <sup>-</sup>	cytochromes <sub>RED</sub>	+290 to +80
MTT + 2H <sup>+</sup> + 2e <sup>-</sup>	MTT <sub>RED</sub>	-110
FMN + 2H <sup>+</sup> + 2e <sup>-</sup>	FMNH <sub>2</sub>	-210
FAD + 2H <sup>+</sup> + 2e <sup>-</sup>	FADH <sub>2</sub>	-220
NAD + 2H <sup>+</sup> + 2e <sup>-</sup>	NADH + H <sup>+</sup>	-320
NADP + 2H <sup>+</sup> + 2e <sup>-</sup>	NADPH + H <sup>+</sup>	-320

A substance on the left side of one of the half reactions can be expected to oxidize a substance on the right side only if the latter has a more negative E<sub>0</sub> than the former. Using this rule, it is evident that MTT will be reduced by FMNH<sub>2</sub>, FADH<sub>2</sub>, NADH, NADPH, but will not be reduced by cytochromes. It is further evident that UptiBlue will be reduced not only by each of these, but by the cytochromes as well.

The importance of this lies in the fact that the flow of electrons may be interrupted by the introduction of an artificial electron-acceptor (redox indicator) with an oxidation-reduction potential intermediate between those of any two members of the electron transport chain. Thus, whenever a substrate is oxidized in the presence of a tetrazolium salt (MTT), the released electrons will not be transported through the usual sequence of cytochromes, but will be trapped. This shuts down the respiratory chain. UptiBlue, on the other hand, is intermediate only between final reduction of O<sub>2</sub> and cytochrome oxidase (Cyt.a<sub>3</sub>). UptiBlue is reduced by the removal of oxygen and its replacement by hydrogen. UptiBlue may substitute for molecular oxygen for any of the oxidoreductases which routinely utilize molecular oxygen as an electron acceptor.

Data may be collected using either fluorescence-based or absorbance-based instrumentation. Fluorescence is monitored at 530-560nm excitation wavelength and 590 nm emission wavelength (Fig.1a). Absorbance is monitored at 570 nm and 600 nm (Figure 1b).

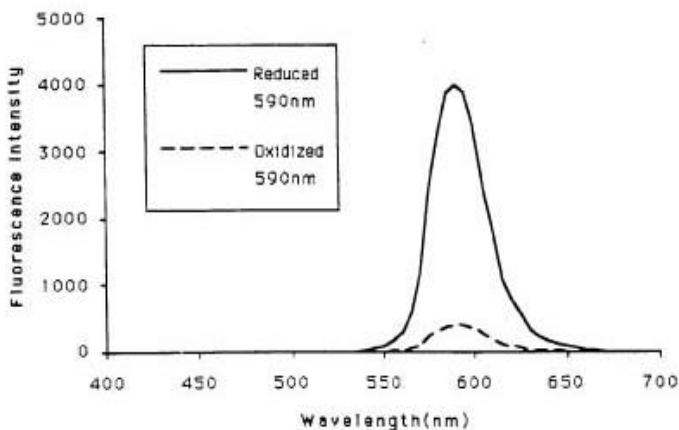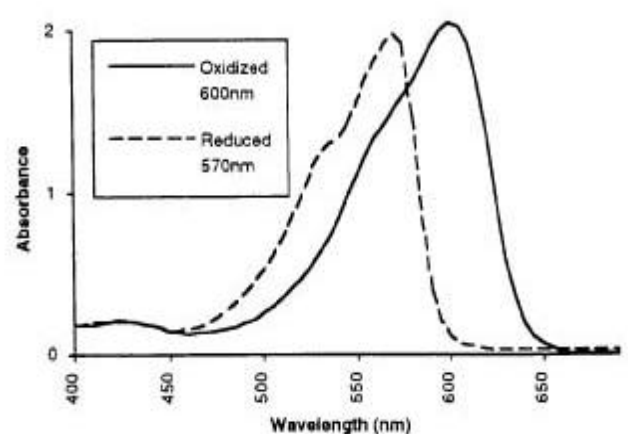
**Figure 1a UptiBlue Fluorescence Emission Spectra**

**Figure 1b UptiBlue Absorbance**

## Advantages

The UptiBlue Assay offers many advantages over conventional cell or radioactively-labelled incorporation assays :

Features	Benefits
Fluorescent/Colorimetric	Allows Choice of detection method
Water soluble	No extraction required
Works on suspension or attached cell lines	No centrifugation required
Fewer steps	Time saving/easily adaptable to automation
Stable	Allows for continuous cell growth monitoring, kinetic studies, incubation time of days
Non-toxic to cells	Less likely to interfere with normal metabolism
Non-toxic to technician	Safe, disposable, less regulation

## Storage Conditions


UptiBlue should be stored in the dark, since the compound is light sensitive (Table 2). The product may be stored for 12 months at room temperature. This expiration date is given on the product label. If shelf life beyond 12 months is desired, storage at 2-8°C increases shelf life to 20 months. UptiBlue may also be frozen at -70°C indefinitely. Because the indicator is a multicomponent solution, it is recommended that frozen UptiBlue be warmed to 37°C and shaken to ensure all components are completely in solution.

### Table 2. UptiBlue Stability

Average absorbance at 600 nm is presented for each month tested. UptiBlue was packaged in amber bottles. Light exposure was continuous at a level of approximately 100 lumens.

Measurements were made on a Perkin Elmer UV/VIS Spectrophotometer Lamda 2 model with a 1 cm path length.

Storage Condition	MONTH									
	0	1	2	3	4	5	6	10	12	
Room Temp Light	0,909	0,925	0,925	0,896	0,901	0,853	0,853	0,820	0,805	
Room Temp Dark	0,909	0,944	0,937	0,944	0,949	0,954	0,954	0,966	0,966	

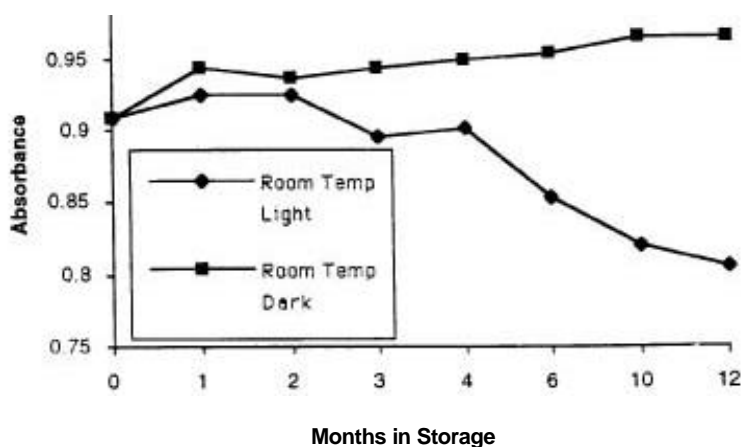

## Quality Control Testing Method for UptiBlue Viable Cell Counting Kit

### Materials and Equipment

- 0,1 M potassium phosphate buffer, pH 7.4
- 10 ml test tube
- pipettor capable of accurately dispensing 0,4 ml
- plate reader with one of the following filters : 540, 570, 600, 630 nm.
- Dynatech flat bottom plate

### Procedure

1. Pipet 0,4 ml of UptiBlue into a test tube.
2. Dilute to 10 mls with phosphate buffer.
3. Mix well.
4. Pipette 100  $\mu$ l into each well of a clear, flatbottom microplate
5. (note : there will enough solution to fill 10 columns in the plate).
6. Read absorbance at appropriate wavelenghts.

- Expected Results :

Wavelength (nm)	Average Absorbance (Standard Deviation)
540	0,145 (0,002)
570	0,225 (0,003)
600	0,313 (0,004)
630	0,116 (0,002)

\* Absorbance Values may be affected by the type of plate (whether round or flat bottom) and the plate manufacturer.

The reduced form (red) of UptiBlue is very unstable in water. For this reason, it is difficult to recommend a standard test for the reduced form. However, the reduced form is very stable in media. To determine the absorbance/fluorescence to be expected from the reduced form (red) for a particular experiment, it is suggested that 1X UptiBlue be made up in the media intended for use in an autoclavable container. Reduce this preparation by autoclaving for 15 minutes. Remove from the autoclave and allow to cool to room temperature. Swirl the solution several times and pipette 100  $\mu$ l into the wells of a flat bottom microtiter plate. Measure the absorbance at the proper wavelengths.

Table 3 a-b presents the absorbance values for the oxidized and reduced forms of the indicator in several commonly used culture media. No QC protocol is recommended for fluorescence since fluorescence units are arbitrary and the scale used varies widely from one instrument to another. From table 3 a-b it should be apparent that the reduced form of UptiBlue is highly fluorescent. When attempting to measure very small changes in reduction, fluorescence measurements will produce greater sensitivity.

### Table 3a-b. Absorbance values for oxidized/reduced forms of UptiBlue for commonly used culture media.

Powdered media was obtained from Sigma and prepared according to their instructions. All media contained phenol red. The required amount of sodium bicarbonate was added to each media, and each was pH adjusted to 7.4 with 1N HCl or 1N NaOH. UptiBlue was added to each media which was then split into 2 samples. One sample of each media was autoclaved for 15 minutes to produce the reduced form. The medias were dispensed into a flat bottom Dynatech plate (100  $\mu$ l per well) and absorbance read at 540, 570, 600, 630 nm on a Cambridge Technologies plate reader. Fluorescence measurements were made on a Cambridge Inc. (Watertown, Mass.) Model 7620 Microplate fluorometer – settings were : bottom reading, light source setting 12, no max AFU, excitation : 530, emission : 560, gain/16.

**a. Absorbance Values**

Powdered Media	Sigma Product #	Wavelength (nm)							
		540		570		600		630	
		Ox	Red	Ox	Red	Ox	Red	Ox	Red
BME EBSS	B9638	0,610	1,207	0,853	1,502	0,845	0,244	0,261	0,177
BME HBS	B9763	0,468	1,087	0,705	1,403	0,817	0,154	0,254	0,097
McCoy's 5A	M4892	0,520	1,133	0,740	1,421	0,756	0,250	0,236	0,183
MEM EBSS	M0268	0,582	1,186	0,819	1,483	0,820	0,235	0,252	0,168
MEM HBSS	M4642	0,480	1,066	0,713	1,383	0,811	0,145	0,251	0,088
Nut Mix F-10	N6635	0,361	0,784	0,583	1,117	0,798	0,138	0,248	0,091
Nut Mix F-12	N6760	0,374	0,796	0,604	1,135	0,822	0,137	0,255	0,085
RPMI 1640	R6504	0,431	0,928	0,659	1,250	0,795	0,161	0,248	0,101

**b. Fluorescence Values**

Powdered Media	Sigma Product #	Fluorescence Units	
		Oxidized	Reduced
BME EBSS	B9638	1926	55676
BME HBSS	B9763	3840	60256
McCoy's 5A	M4892	2640	50545
MEM EBSS	M0268	2377	54493
MEM HBSS	M4642	4194	59202
Nut Mix F-10	N6635	2472	70092
Nut Mix F-12	N6760	5232	68132
RPMI 1640	R6504	6472	58796

**General Procedure : Determining Length of Incubation Time and Plating density for a Cell Line**

The two variables which most affect the response of cells to UptiBlue are length of incubation time and number of cells plated. It is recommended that the plating density and incubation time be determined for each cell line using the following procedure :

1. Harvest cells which are in log phase growth stage and determine cell count. Plate cells at various densities, some dilutions being above and below the cell density expected to be used.
2. Aseptically add UptiBlue in an amount equal to 10% of the culture volume.
3. Return cultures to incubator. Remove the plate and measure fluorescence/absorbance each hour following plating for the first 6-8 hours. It is also recommended that the plate remain in incubation overnight and measurements be made the following day at 24 hours. Two kinds of information can be obtained from this data : (1) for any given incubation time selected, the range in cell density can be determined for which there is a linear response relating cell numbers to UptiBlue reduction, and (2) for any given cell density selected, the maximum incubation time can be determined in which the control cells turn the indicator from the oxidized (blue) form to the fully reduced (red) form.
4. Measure absorbance at a wavelength of 570 nm and 600 nm. Or, measure fluorescence with excitation wavelength at 530-560 nm and emission wavelength at 590 nm.
5. To generate the graph for (1), plot the log of cell density on the x-axis and reduction of UptiBlue from absorbance (using equation 5 to be discussed in calculations for absorbance) or fluorescence on the y-axis. (fig.2) To generate the graph for (2), plot the number of hours incubated on the x-axis and reduction of UptiBlue on the y-axis (Fig.3) A sample data set and the resulting graphs are presented below.

NT-UP66941  
Sample Data Set

Time	Blue in Media	a549 cells/ml				p388 cells/ml			
		10 000	5 000	1 000	500	10 000	5 000	1 000	500
T0	0,336	0,372	0,348	0,338	0,336	0,328	0,332	0,334	0,342
T4	0,334	0,540	0,432	0,352	0,339	0,335	0,335	0,333	0,340
T4	0,333	0,590	0,489	0,365	0,346	0,346	0,339	0,332	0,339
T5,5	0,321	0,573	0,511	0,366	0,344	0,344	0,335	0,325	0,331
T20	0,322	0,486	0,518	0,510	0,438	0,434	0,381	0,328	0,332

Absorbance values at 570 nm after blanking with media only. a549 is a monolayer culture, p388 is a suspension cell line.

Time	Blue in Media	a549 cells/ml				p388 cells/ml			
		10000	5000	1000	500	10000	5000	1000	500
T0	0,441	0,496	0,459	0,443	0,439	0,422	0,433	0,438	0,451
T4	0,440	0,267	0,349	0,421	0,425	0,404	0,422	0,435	0,448
T4	0,432	0,162	0,265	0,397	0,411	0,391	0,414	0,432	0,444
T5,5	0,424	0,135	0,211	0,377	0,397	0,377	0,404	0,424	0,435
T20	0,412	0,112	0,102	0,180	0,271	0,253	0,337	0,415	0,423

Absorbance values at 600 nm after blanking with media only. a549 is a monolayer culture, p388 is a suspension cell line.

Using absorbance data from the sample data set, percent reduction of UptiBlue was calculated using equation 5 :

Time	a549 cells/ml				p388 cells/ml			
	10 000	5 000	1 000	500	10 000	50 00	1 000	500
T0	5,7	6,0	6,1	6,3	7,0	6,3	6,0	5,9
T2	65,7	35,4	11,5	8,6	10,6	8,3	6,3	5,9
T4	89,9	57,6	17,3	11,9	14,5	10,2	6,6	6,3
T5,5	91,8	70,0	20,4	13,6	16,2	10,9	6,4	6,1
T20	80,7	88,3	76,2	49,6	51,3	29,5	8,4	8,1

These values are plotted to produce figures 2-3

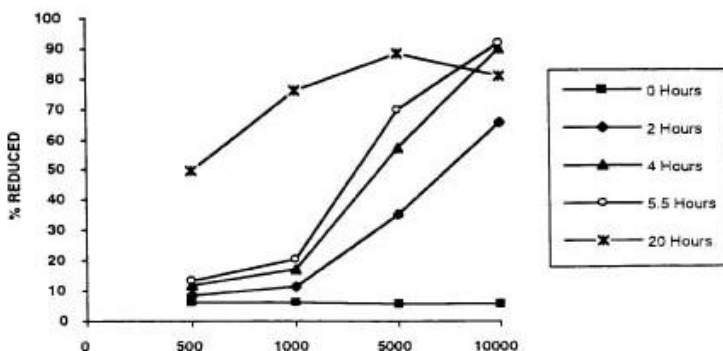


Fig. 2A. A549

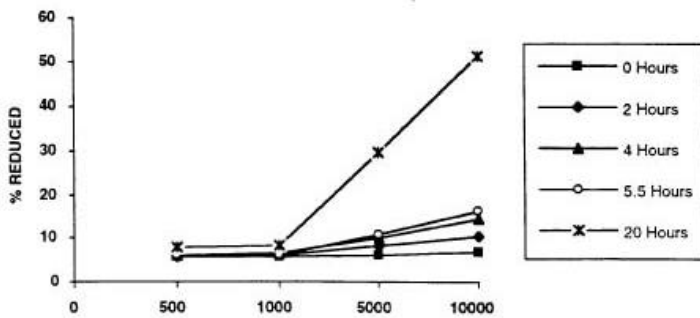


Fig. 2B. P388

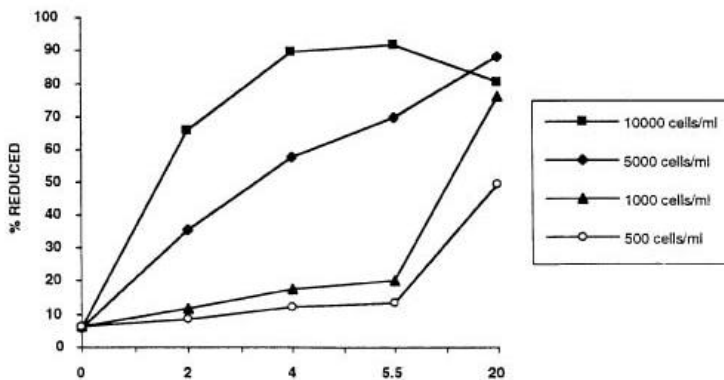


Fig. 3A. A549

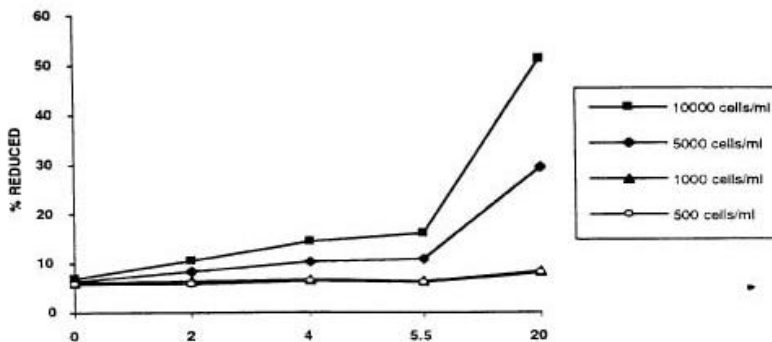


Fig. 3B. P388

From fig. 2A, if, for example, the desired incubation time with UptiBlue were 4 hours, any plating density from 500 to 10,000 cells/ml could be used and expected to produce a reaction within the linear range for UptiBlue for that incubation period. However, if the intent is to incubate for 20 hours, the reaction could only be expected to be within the linear range if plated at 500-1000 cells/ml for this cell line.

On the other hand, with P388, even with an initial plating density up to 10,000 cells per well, data is within range and UptiBlue has only been reduced by 50% at that point. This indicates cells could be incubated with UptiBlue, even when plated at this high density for up to 2 days. If a shortened exposure is desired, then the initial plating density should be increased. If the goal is to continue the experiment for more than 2 days, the initial plating density should be decreased. Similar information is gained from examination of fig. 3a-b.

**Note**

With high cell numbers or extended incubation time (days), you will reach a point where the red form stops increasing and begins to decline. The absorbance/fluorescence level drops with a corresponding clearing of the red colour. This is demonstrated in Fig. 3 A when plated at 10,000 cells/ml and incubated for more than 4 hours.

Microbial contaminants will also reduce UptiBlue and will yield erroneous results if contaminated cultures are tested by this method.

Samples with protein concentrations equivalent to 10% foetal bovine serum do not interfere with the assay (8). Further, there is no interference from the presence of phenol red in the growth medium. The presence of phenol red merely shifts the values approximately 0,03 units higher (see table 4).

**Table 4. Effect of phenol red on absorbance values at 570 nm.**

Absorbance value for various levels reduction of UptiBlue in RPMI 1640 w/MOPS with and without phenol red, pH 7.0. 100 µl per well, Dynatech flat bottom plate.

	Media	Blue	% REDUCED				Red
			10	30	60	90	
RPMI 1640	0,032	0,47	0,52	0,61	0,73	0,85	0,88
RPMI 1640 w/phenol red	0,061	0,53	0,54	0,64	0,76	0,88	0,91

**Effect of Storage of Plates on Measurements**

Many investigators find that they may be able not to read plates on the day an experiment is performed. It is recommended that plates be refrigerated and read within 1-3 days. Plates can be wrapped in foil or plastic wrap to prevent evaporation. Table 5 presents absorbance data for the oxidized (blue) and reduced (red) forms of UptiBlue for plates which were read on day 1, stored overnight refrigerated and read again on days 2 and 3. Data is presented for some different plate types and plate manufacturers to illustrate the effect these variables have on absorbance values obtained.

Table 6 gives the fluorescence values for the same plate. If plates are stored refrigerated and fluorescence measurements are being used, keep in mind that fluorescence measurements are influenced by temperature (see table 7). If measurements are normally taken at 37°C, then plates should be warmed to that temperature before reading.

**Table 5. Effect of Storage Plates on Absorbance**

100 µl of RPMI 1640 w/MOPS 7.0 no phenol red.

		ABSORBANCE							
		BLUE (Oxidized)				RED (Reduced)			
		540 nm	570 nm	600 nm	630 nm	540 nm	570 nm	600 nm	630 nm
Dynatech Flat Bottom	Day 1	(.003)	(.005)	(.007)	(.003)	(.020)	(.027)	(.002)	(0.0)
		.298	.496	.708	.236	.693	1.017	.126	.075
	Day 2	(.003)	(.004)	(.006)	(.002)	(.020)	(.027)	(.008)	(.009)
		.294	.484	.692	.227	.697	1.018	.164	.118
	Day 3	(.003)	(.006)	(.008)	(.003)	(.010)	(.024)	(.008)	(.009)
		.296	.486	.691	.231	.734	1.038	.199	.149



NT-UP66941

Corning Flat Bottom	Day 1	(.002) .210	(.004) .335	(.005) .474	(.002) .169		(.020) .530	(.024) .772	(.003) .137	(.004) .105
	Day 2	(.002) .210	(.003) .329	(.004) .458	(.001) .161		(.020) .580	(.027) .822	(.004) .193	(.004) .159
	Day 3	(.002) .200	(.003) .322	(.005) .444	(.002) .160		(.020) .600	(.035) .823	(.004) .210	(.003) .172
Corning Round Bottom	Day 1	(.001) .380	(.002) .635	(.002) .913	(.001) .300		(.014) .870	(.018) 1.266	(.002) .151	(0,0) .084
	Day 2	(.002) .390	(.003) .641	(.004) .914	(.002) .295		(.011) .850	(.014) 1.241	(.002) .146	(.002) .083
	Day 3	(.004) .390	(.006) .646	(.008) .916	(.004) .302		(.017) .860	(.021) 1.237	(.007) .159	(.006) .094

**Table 6. Effect of Storage of Plates on Fluorescence**

Mesasurements made at room temperature (22°C)

	Fluorescence Units						
	BLUE (Oxidized)				RED (Reduced)		
	Day 1	Day 2	Day 3		Day 1	Day 2	Day 3
Dynatech Flat Bottom	(7.4) 188	(49.6) 265	(42.8) 277		(123) 5821	(127) 5973	(157) 5867
Corning Flat Bottom	(4) 160	(6) 184	(10) 206		(101) 5065	(101) 5124	(130) 4904
Corning Round Bottom	(4.1) 185	(16) 250	(43) 289		(61) 6048	(219) 5904	(110) 5481

Standard Deviation in parentheses (calculated for n = 8)

**Table 7. Effect of Temperature on Fluorescence**

RPMI 1640 w/MOPS pH 7.0 no phenol red, 100 µl per well, Dynatech flat bottom.

Temp.	Fluorescence
37° C	(171,1) 5216
22° C	(157,6) 5867
4° C	(192,1) 6881

## General Discussion of calculations when Using Absorbance.

It is clear from Fig.1b that there is considerable overlap in the absorption spectra of the oxidized and reduced forms of UptiBlue . When there is no region in which just one component absorbs, it is still possible to determine the two substances by making measurements at two wavelengths. The two components must have different powers of light absorption at some points in the spectrum. Since absorbance is directly proportional to the product of the molar extinction coefficient and concentration, a pair of simultaneous equations may be obtained from which the two unknown concentrations may be determined :

$$1) C_{RED} (\epsilon_{RED}) \lambda_1 + C_{OX} (\epsilon_{OX}) \lambda_1 = A \lambda_1$$

$$2) C_{RED} (\epsilon_{RED}) \lambda_2 + C_{OX} (\epsilon_{OX}) \lambda_2 = A \lambda_2$$

To solve for the concentration of each component :

$$3) C_{RED} = \frac{(\epsilon_{OX}) \lambda_2 A \lambda_1 - (\epsilon_{OX}) \lambda_1 A \lambda_2}{(\epsilon_{RED}) \lambda_1 (\epsilon_{OX}) \lambda_2 - (\epsilon_{OX}) \lambda_1 (\epsilon_{RED}) \lambda_2}$$

$$4) C_{OX} = \frac{(\epsilon_{RED}) \lambda_1 A \lambda_2 - (\epsilon_{RED}) \lambda_2 A \lambda_1}{(\epsilon_{RED}) \lambda_1 (\epsilon_{OX}) \lambda_2 - (\epsilon_{OX}) \lambda_1 (\epsilon_{RED}) \lambda_2} = 100$$

To determine the percent reduction of UptiBlue :

$$\% \text{ Reduced} = \frac{C_{RED \text{ Test Well}}}{C_{OX \text{ Negative Control Well}}}$$

$$5) = \frac{(\epsilon_{OX}) \lambda_2 A \lambda_1 - (\epsilon_{OX}) \lambda_1 A \lambda_2}{(\epsilon_{RED}) \lambda_1 A' \lambda_2 - (\epsilon_{RED}) \lambda_2 A' \lambda_1} \times 100$$

To calculate the percent difference in reduction between treated and control cells in cytotoxicity/proliferation assays :

$$6) \frac{(\epsilon_{OX}) \lambda_2 A \lambda_1 - (\epsilon_{OX}) \lambda_1 A \lambda_2 \text{ of test agent dilution}}{(\epsilon_{OX}) \lambda_2 A^\circ \lambda_1 - (\epsilon_{OX}) \lambda_1 A^\circ \lambda_2 \text{ of untreated positive growth control}} \times 100$$

Where

$C_{RED}$  = concentration of reduced form UptiBlue (RED)

$C_{OX}$  = oxidized form of UptiBlue (BLUE)

A = absorbance of test wells

A' = absorbance of negative control well. The negative control well should contain media + UptiBlue but no cells.

A° = absorbance of positive growth control well.

$\lambda_1$  = 570 nm (540 nm may also be used).

$\lambda_2$  = 600 nm (630 may also be used).

**The key** equations to use are (5) and (6), equation 5 for calculating the percent reduction from the blue oxidized form and equation 6 for calculating percent difference between treated and control cells in cytotoxicity/proliferation experiments. Equations 1-4 are included only for completeness of the discussion.

Blanking of the plate reader should be done with a well containing media only. Table 3 contains the necessary values for solving the equation stated above. Table 4 lists typical absorbance values for different percent reduction of UptiBlue.

**TABLE 3. Molar Extinction Coefficients for UptiBlue .**

Wavelength (λ)	ε <sub>RED</sub>	ε <sub>Ox</sub>
540 nm	104,395	47,619
570 nm	155,677	80,856
600 nm	14,652	117,216
630 nm	5,494	34,798

**TABLE 4. Typical Absorbance Values for Microplate Reader.** Test wells contain 100 µl RPMI 1640 w/MOPS with no phenol red at pH 7.0 with UptiBlue, containing the reduced form present in known amounts.

% Reduced UptiBlue	Wavelength (nm)			
	540	570	600	630
0 Fully Oxidized (Blue)	0,29	0,47	0,67	0,22
10	0,29	0,52	0,62	0,21
20	0,35	0,56	0,57	0,19
30	0,38	0,61	0,50	0,18
40	0,41	0,64	0,45	0,16
50	0,45	0,69	0,39	0,14
60	0,48	0,73	0,33	0,13
70	0,51	0,77	0,28	0,11
80	0,54	0,81	0,23	0,10
90	0,57	0,85	0,16	0,08
100 Fully Reduced (Red)	0,60	0,88	0,10	0,06
Media	0.03	0.03	0.03	0.03

**Example Calculation 1 :**

$$\lambda_1 = 570, \lambda_2 = 600$$

$$(\epsilon_{Ox})\lambda_2 = 117,216$$

$$(\epsilon_{RED})\lambda_2 = 80,856$$

$$(\epsilon_{RED})\lambda_1 = 155,677$$

$$(\epsilon_{RED})\lambda_2 = 14,652$$

$$A\lambda_1 = 0,61$$

Observed absorbance reading for test well

$$A\lambda_2 = 0,42$$

Observed absorbance reading for test well

$$A'\lambda_2 = 0,64$$

Observed absorbance reading for negative control well

$$A'\lambda_1 = 0,44$$

Observed absorbance reading for negative control well

To calculate percent reduced using equation 5 :

Percent reduced

$$= \frac{(\epsilon_{Ox})\lambda_2 A\lambda_1 - (\epsilon_{Ox})\lambda_1 A\lambda_2}{(\epsilon_{RED})\lambda_1 A'\lambda_2 - (\epsilon_{RED})\lambda_2 A'\lambda_1} \times 100$$

$$= \frac{(117,216)(,61) - (80,586)(,42)}{(155,677)(,64) - (14,652)(,44)} \times 100$$

$$= \frac{71,502 - 33,846}{99,633 - 6,446} = \frac{37,655}{99,186} = .404 \times 100 = 40 \%$$

**Example Calculation 2 :**

$\lambda_1 = 570, \lambda_2 = 600$   
 $(\epsilon_{0x}) \lambda_2 = 117,216$   
 $(\epsilon_{0x}) \lambda_1 = 80,856$   
 $(\epsilon_{RED}) \lambda_1 = 155,677$   
 $(\epsilon_{RED}) \lambda_2 = 14,652$   
 $A \lambda_1 = 0,65$  Observed absorbance reading for test well  
 $A \lambda_2 = 0,36$  Observed absorbance reading for test well  
 $A^o \lambda_2 = 0,78$  Observed absorbance reading for positive growth control  
 Well  
 $A^o \lambda_1 = 0,19$  Observed absorbance reading for positive growth control  
 Well

To calculate the percent difference in reduction between treated and control cells in cytotoxicity/proliferation assays using equation 6 :

Percent difference in reduction

$$= \frac{(\epsilon_{0x})\lambda_2 A\lambda_1 - (\epsilon_{0x})\lambda_1 A\lambda_2 \text{ of test agent dilution}}{(\epsilon_{0x})\lambda_2 A^o\lambda_2 - (\epsilon_{0x})\lambda_1 A^o\lambda_1 \text{ of untreated positive growth control}} \times 100$$

$$= \frac{(117,216)(.65) - (80,856)(.36)}{(117,216)(.78) - (80,856)(.19)} \times 100$$

$$= \frac{76,190 - 29,010}{91,428 - 15,311} \times 100$$

$$= \frac{47,199}{76,116} \times 100 = 62 \%$$

This would indicate that the amount of reduction in the test well is only 62% of that in the control well, or put another way, that growth in the test well is inhibited by 38% when compared to that of the control.

**UptiBlue Reduction Curves**

Examples of reduction curves are included to demonstrate the usefulness of the UptiBlue Assay for measuring cell proliferation. (Fig. 4)

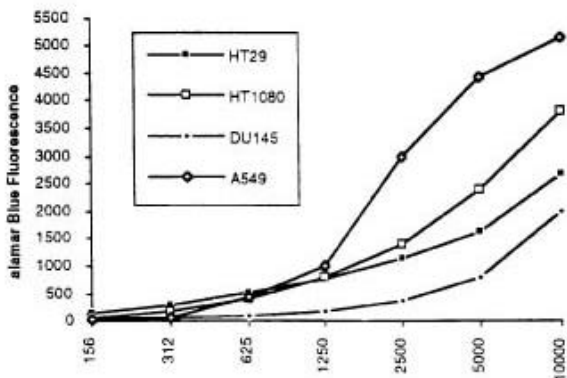


Figure 4 Detection of cell Growth of Four Cell Lines using UptiBlue

NT-UP66941

UptiBlue is especially well suited for kinetic studies that involve monitoring cell growth for extended exposure periods. The stability of UptiBlue allows the investigator to add the indicator at the beginning of the experiment and continue to follow reduction by the cells for several days. (Figure 5)

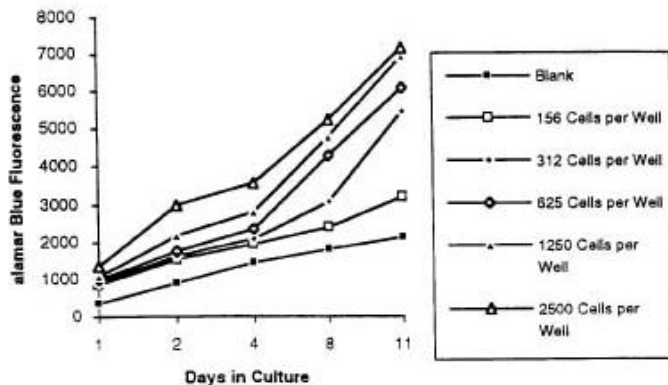


Figure 5 Continuous Incubation of Cell Line A549 with UptiBlue

## Example Procedure of Cytotoxicity assay

### Preparation of cells for Testing

1. Harvest an appropriate cell line by trypsinization and subsequent trypsin inhibitor treatment.
2. Centrifuge cells, resuspend in growth medium and count.
3. Calculate the total cell number and adjust to  $1 \times 10^4$  cell/ml. This is a suggested cell density which has worked in our studies with cancer cell lines. Refer to the article by Alley et al. (3) for suggestions on plating densities and growth medium when working with cancer cell lines and chemotherapeutic agents.
4. Add 250  $\mu$ l of cell suspension to each well. Incubate at 37° in 5% CO<sub>2</sub> atmosphere for the number of days required for the particular cell line to be in log phase (usually 3 days).

### Exposing Cells to Test Agents

1. Prepare appropriate dilutions of test agent in growth media.
2. Aspirate spent growth medium from the wells and add 250 $\mu$ l of each dilution of test agent to the wells.
3. Cover, then return to the incubator for 2 days.
4. After incubation, add 25  $\mu$ l of the indicator to each well. Incubate panels for an additional 3 hours. Panels may then be read spectrophotometrically (absorbance at 570 nm and 600 nm) or spectrofluorometrically (excitation : 530-560 nm ; emission : 590 nm).

### Data Analysis

#### Fluorescence :

1. Calculate percent of untreated control with the following formula :

$$\frac{\text{Fl 590 of test agent dilution}}{\text{Fl 590 of untreated control}} \times 100$$

2. Use semi-log graph paper and plot the percent of untreated control for each dilution of a given test agent on the y-axis vs. the concentration of the test agent on the x-axis.

NT-UP66941


- Determine the LD<sub>50</sub> endpoint from the graph by reading from where the 50 percent point intercepts the Dose Response Curve to the concentration along the x-axis. That concentration is the LD<sub>50</sub> value. (Figure 6)

**Absorbance :**

- Calculate percent of untreated control with the following formula :

$$= \frac{(\epsilon_{0x})\lambda_2 A\lambda_1 - (\epsilon_{0x})\lambda_1 A\lambda_2 \text{ of test agent dilution}}{(\epsilon_{0x})\lambda_2 A\lambda_1 - (\epsilon_{0x})\lambda_1 A\lambda_2 \text{ of untreated positive growth control}} \times 100$$

- Use semi-log graph paper and plot the percent of untreated control for each dilution of a given test agent on the y-axis vs. the concentration of test agent on the x-axis
- Determine the LD<sub>50</sub> endpoint from the graph by reading from where the 50 percent point intercepts the Dose Response Curve to the concentration along the x-axis. That concentration is the LD<sub>50</sub> value. (Figure 6)

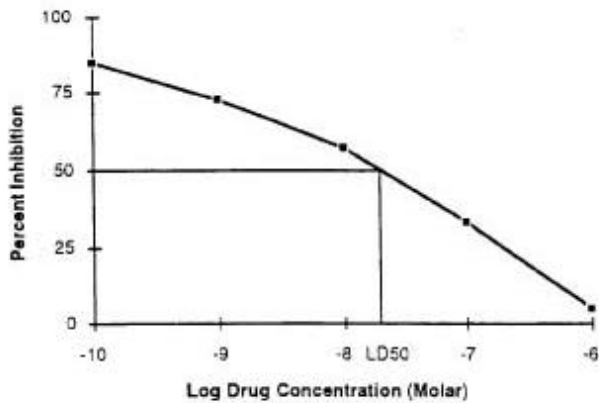


**Figure 6 Determination of Doxorubicin LD<sub>50</sub> using UptiBlue**


**Applications**

Current users of MTT, XTT or neutral red uptake in proliferation/cytotoxicity assays can substitute UptiBlue for each of these tests. This substitution can be done at the time point when you would normally add MTT, XTT or neutral red.


Figures 7-9 are examples from comparison tests with each of these methods.

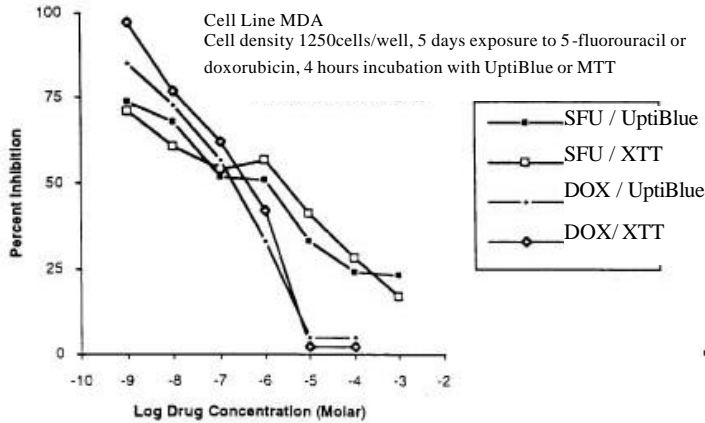

**Figure 7 Comparison of LD<sub>50</sub> Using UptiBlue and MTT**

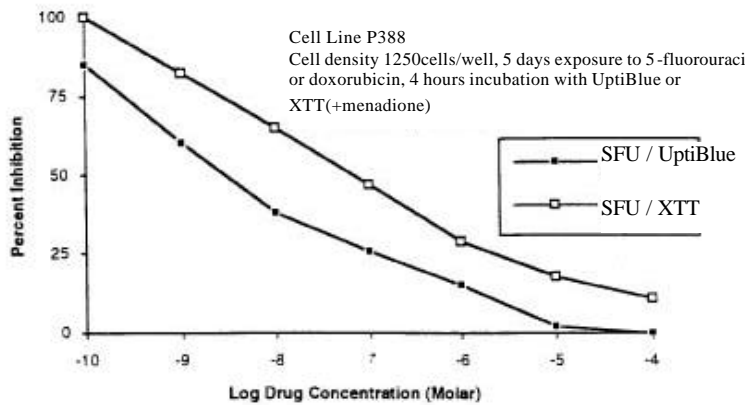

**Figure 8. Comparison of LD<sub>50</sub> Using UptiBlue and XTT LD<sub>50</sub> for UptiBlue is derived from fluorescence measurements**

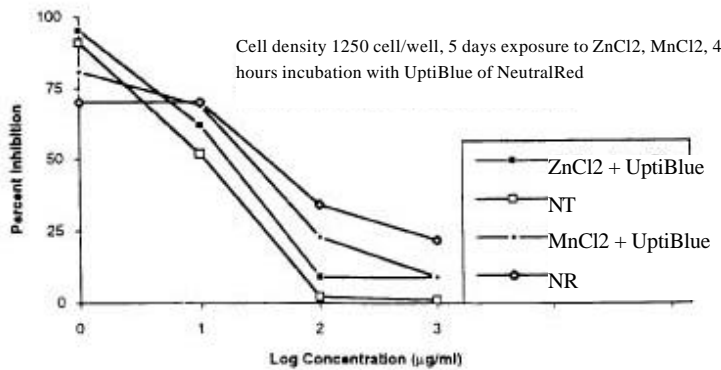

**Figure 9. Comparison of LD<sub>50</sub> using UptiBlue and Neutral Red for Abrasives Applied to Human Epidermal Keratinocytes LD<sub>50</sub> for UptiBlue is derived from fluorescence measurements.**