

KOBRA[®] CELL Art. No. K01

Instruction Manual



● Contents

1.	Introduction to the KOBRA® CELL	4
1.1	Intended Use	4
1.2	Derivatisation Reaction	5
2.	KOBRA® CELL Contents.....	6
2.1	Contents.....	6
2.2	Reagents Required	6
2.3	Accessory Products	6
2.4	Hazards	6
2.5	Decontamination.....	6
3.	Assembly of the KOBRA® CELL.....	7
3.1	Components	7
3.2	Precautions When Assembling.....	8
3.3	Assembly	8
4.	KOBRA® CELL Installation.....	9
4.1	Power Pack Notes	9
4.2	Power Pack Technical Specification	9
4.3	Installation Guide	10
4.4	Installation Instructions	11
4.5	Precautions to Take Before Running the KOBRA® CELL.....	12
4.6	Operating the KOBRA® CELL	12
4.7	Error LED	12
5.	HPLC Conditions	13
6.	KOBRA® CELL Maintenance	14
6.1	Storage.....	14
6.2	Daily Cleaning	15
6.3	Monitoring Performance.....	16
6.4	Changing the Membrane.....	16

7.	Troubleshooting.....	17
7.1	No Peaks Observed.....	17
7.2	Reduced / No Aflatoxin B1 and G1 Peaks Observed	18
7.3	Abnormal Peaks Observed.....	19
7.4	Error Light Observed on Power Pack	20
7.5	Leakage from KOBRA® CELL	20
7.6	KOBRA® CELL Causing High Pressure on the HPLC	20
7.7	Storage Plug Broken Off or Lodged in KOBRA® CELL	20
8.	General Information	21
8.1	Quality	21
8.2	Technical Support.....	21
8.3	Warranty.....	21

1 Introduction to KOBRA[®] CELL

1.1 Intended Use

Confirmation of the presence of aflatoxins in a sample by HPLC requires derivatisation of aflatoxins B1 and G1 in order to enhance their natural fluorescence and make them more easily detected. Previously, the only options available for derivatising aflatoxins involved the use of trifluoroacetic acid (TFA), pyridinium bromide perbromide (PBPB) or iodine. All of these methods have significant limitations which can be overcome by use of the KOBRA[®] CELL.

Pre-column derivatisation with TFA requires the solution containing aflatoxin to be blown to dryness under a stream of nitrogen, potentially leading to a loss of toxin. Further limitations are that the reaction takes 30 minutes at 50 °C, and the TFA reagent is itself corrosive and harmful to handle.

The post-column PBPB method involves addition of the diluted reagent into the eluate from the HPLC column. The limitations of this method are the use of a second pump and the difficulty in dissolving the PBPB as well as the hazardous nature of the reagent.



Post-column derivatisation with iodine also has some limitations including the need for a second pump, water bath or oven which can be expensive. It is necessary to clean the equipment regularly in order to avoid iodine crystals forming inside the reaction coil. Finally, the iodine must be prepared fresh each day due to its unstable nature.

The KOBRA[®] CELL overcomes the problems relating to alternative derivatisation procedures. It is an electrochemical cell connected to an HPLC system downstream from the HPLC column and in line with the column effluent and the fluorescence detector. The KOBRA[®] CELL generates a reactive form of bromine for derivatisation of aflatoxins B1 and G1, resulting in enhanced fluorescence and thus more sensitive detection.

The KOBRA[®] CELL is used by hundreds of labs around the world and is mentioned in several EU and other international standard methods.

1.2 Derivatisation Reaction

The aflatoxins and the mobile phase enter the KOBRA® CELL and the electrochemical reaction occurs generating the reactive form of bromine. The reaction between the reactive bromine and the aflatoxins must take place before the derivatised aflatoxins enter the fluorescence detector. Hence the length of the reaction coil is critical. A minimum reaction time of 4 seconds is required.

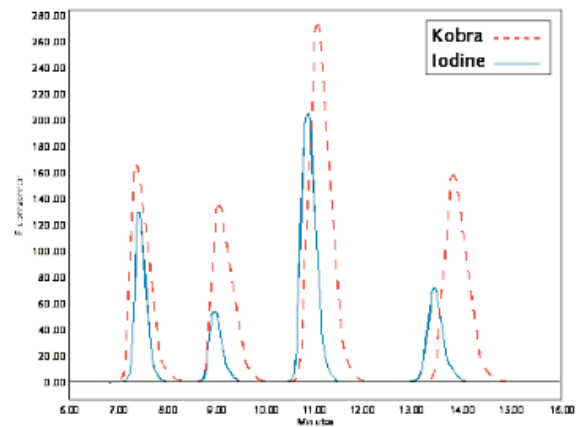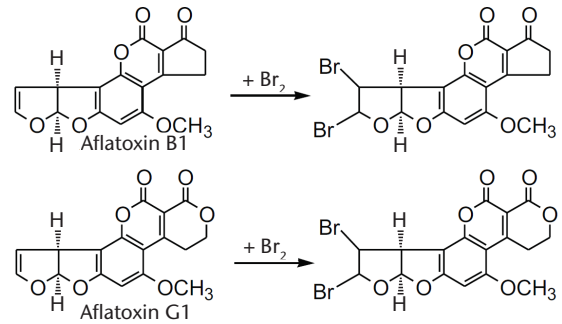

2 KOBRA® CELL Contents

2.1 Contents

- 1 x KOBRA® CELL
- 1 x Power Pack (including 1 red and 1 black connection lead)
- 1 x Electrical adapter (with various adapters)
- 1 x 1 m length of 0.5 mm ID PEEK™ tubing
- 1 x Spare membrane

2.2 Reagents Required

- Distilled Water
- Methanol
- Potassium Bromide
- Nitric Acid

2.3 Accessory Products

- 1 x KOBRA® CELL Installation Pack (Product Code: K03)*
- Guard column: Inertsil ODS-3, 5 µm, 4 mm x 10 mm or equivalent

* Available from R-Biopharm. Please contact your local distributor for further information.

2.4 Hazards

Mycotoxins are very hazardous substances. Only laboratories equipped to handle toxic materials and solvents should perform analyses. Suitable protective clothing, including gloves, safety glasses and lab coats should be worn throughout the analysis.

2.5 Decontamination

Prior to disposal, excess standard solutions should be treated with at least one-tenth their volume of 5 % sodium hypochlorite. Labware and contaminated waste should be immersed in 5 % sodium hypochlorite solution for 30 minutes followed by the addition of 5 % acetone for 30 minutes. Flush with copious amounts of water before disposal. After decontamination labware should be thoroughly washed. Incinerate waste if regulations permit.

3 Assembly of the KOBRA[®] CELL

3.1 Components

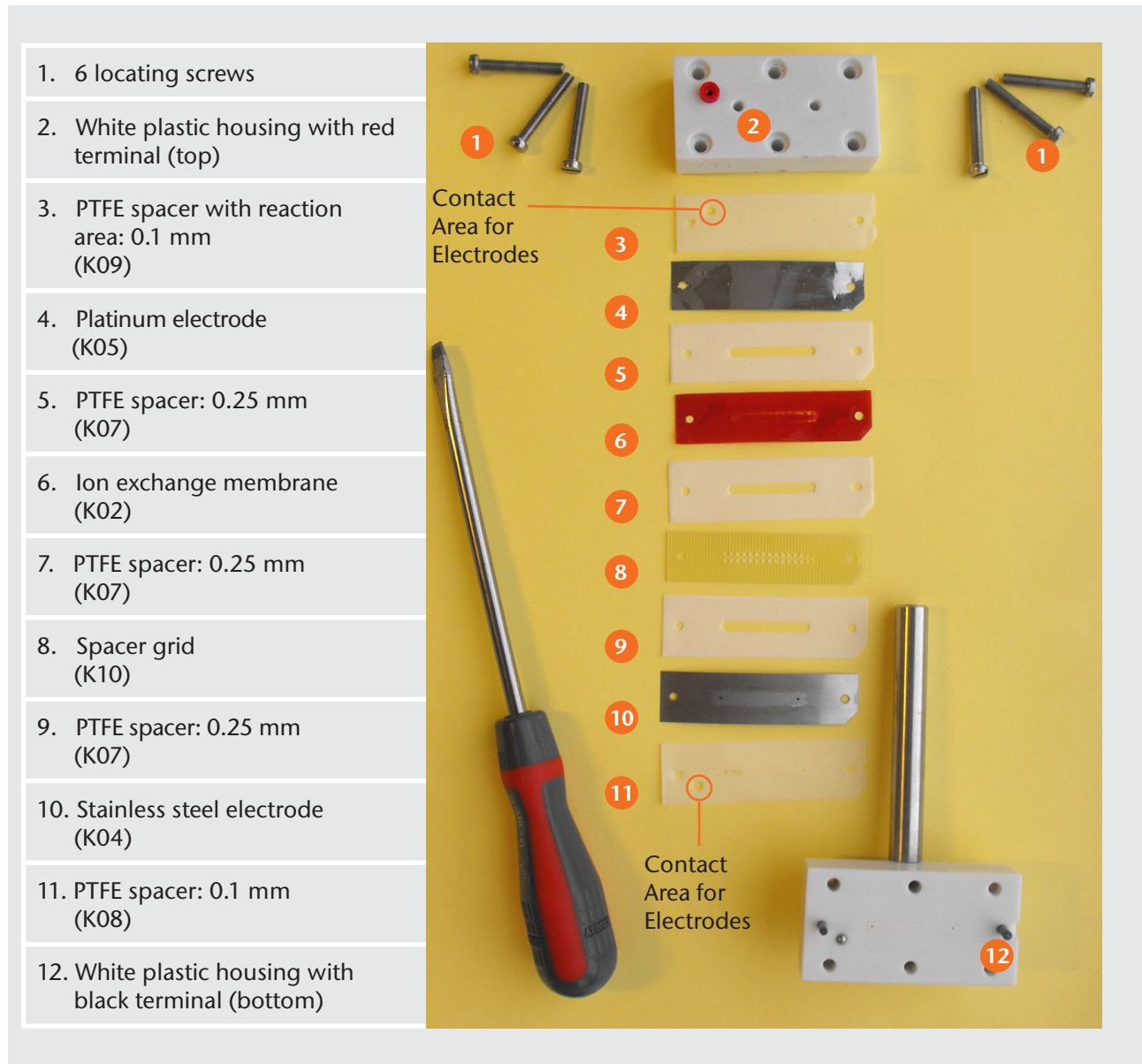


Fig 1. KOBRA[®] CELL Components

3.2 Precautions When Assembling

- To reassemble the KOBRA® CELL it is important to tighten the locating screws at opposite corners before tightening the middle screws. This will ensure that the KOBRA® CELL is correctly assembled and will not leak.
- Do not over tighten locating screws as this will damage the threads.
- Never use metal nuts or ferrules directly with the KOBRA® CELL.
- Use a guillotine tubing cutter to prevent distortion of the diameter of the PEEK™ tubing used when assembling the KOBRA® CELL.

3.3 Assembly

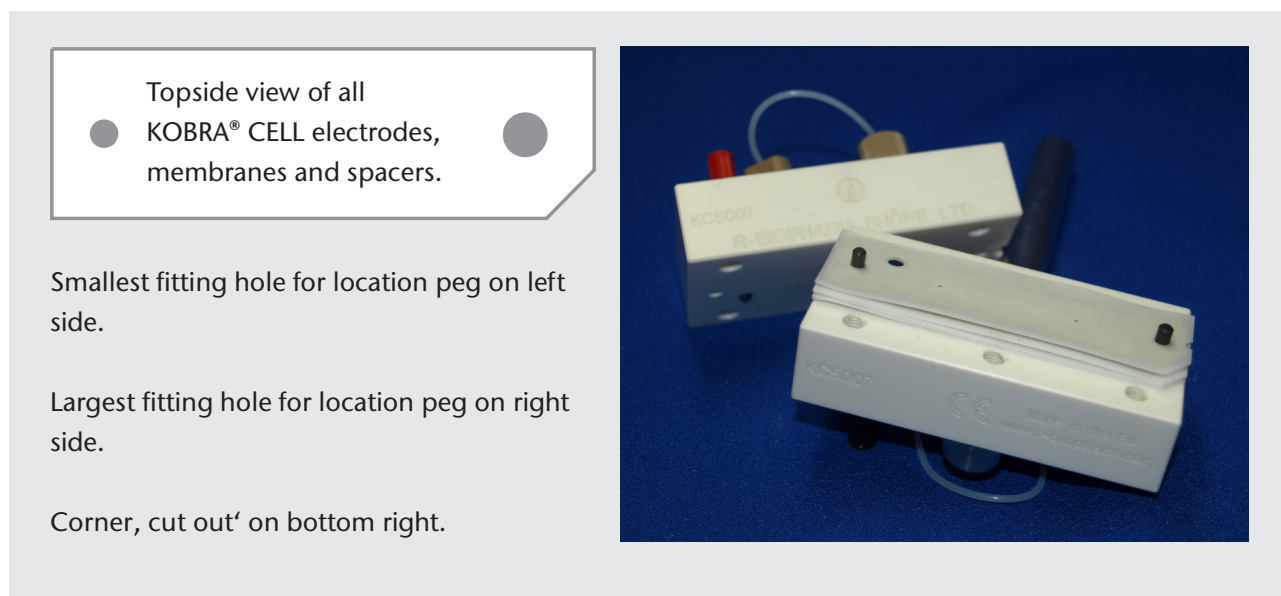


Fig 2. Assembly of KOBRA® CELL

4 KOBRA® CELL Installation

4.1 Power Pack Notes

It should be noted that the KOBRA® CELL power pack has been improved and as such is of a slightly different design. The new power pack is still 9 V and is supplied with an adapter for use in the UK, USA, Europe and Australia. The plug can be converted by depressing the button located at the top of the adapter using a screwdriver. This releases the adapter that is in position. Slide to remove. To attach the required adapter place on the two pins and slide to connect. The power pack can now only be used at 100 μ A and there is a light to show it has been switched on (top right). There is also an error light to indicate if there is a problem (bottom left).

4.2 Power Pack Technical Specification

Power Supply	Primary 110 - 240 V Secondary 9 V
Exit Power Current Source	100 μ A
Maximum Output Power	9 V
Loading Time NICAD 9 V	Approximately 1 hour 15 minutes.
Accuracy	0.5 % (This is an improvement on the previous model which was 2 %.)
Working Time with Battery	20 hours
Approval	CE approval conforms to European Directives 73/23/EEC and 89/336/EEC.

4.3 Installation Guide

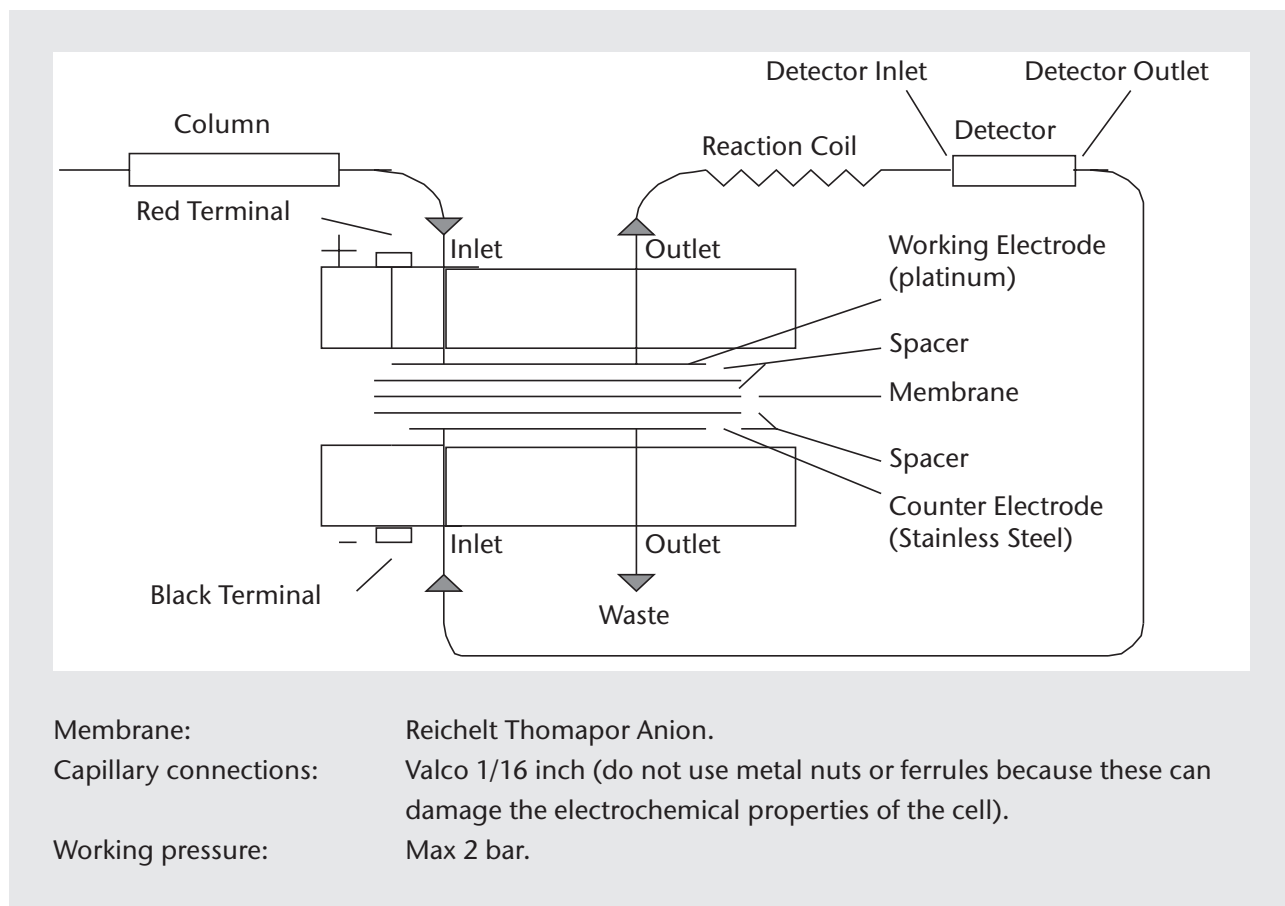


Fig 3. KOBRA® CELL Installation

4.4 Installation Instructions

Caution: Do not over tighten. Cell connections should be made with plastic ferrules, not metal.

1. Unpack the KOBRA® CELL and check that all of the components are present.
2. Add 119 mg of potassium bromide and 350 µl of 4 M nitric acid to 1 l of mobile phase.
3. Immerse the tubing from the HPLC pump into the mobile phase.
4. Disconnect the storage plugs and capillary tubing from the KOBRA® CELL. Connect the HPLC column discharge to the KOBRA® CELL inlet using plastic ferrules (hand tightened) and a length of PEEK™ tubing supplied with the kit.
5. Connect the KOBRA® CELL outlet using ferrules (hand tightened) and a length of PEEK™ tubing to the detector inlet according to the table below.
6. Connect the detector outlet to the KOBRA® CELL inlet using plastic ferrules (hand tightened) and PEEK™ tubing.
7. Connect the KOBRA® CELL outlet to waste using plastic ferrules (hand tightened) and PEEK™ tubing.
8. Switch on the HPLC pump and flush mobile phase through the system for approximately 5 minutes.
9. Connect the KOBRA® CELL power pack to the cell, red lead to red terminal and black lead to black terminal.
10. After the HPLC has been allowed to stabilise (approximately 30 minutes) the KOBRA® CELL is ready to use.

Caution: Do not switch on without mobile phase flowing through the KOBRA® CELL or the membrane will be damaged.

Caution: Use a guillotine tubing cutter to prevent distortion of the diameter of the tubing.

Determination of Reaction Coil / PEEK™ Tubing Length							
Flow Rate		0.5 ml/min	0.6 ml/min	0.7 ml/min	0.8 ml/min	0.9 ml/min	1.0 ml/min
Internal diameter of reaction coil / PEEK™ tubing (id)	0.2 mm	106.1 cm	127.3 cm	148.5 cm	169.8 cm	191.0 cm	212.2 cm
	0.25 mm	67.9 cm	81.5 cm	95.1 cm	108.6 cm	122.2 cm	135.8 cm
	0.4 mm	26.5 cm	31.8 cm	37.1 cm	42.4 cm	47.7 cm	53.1 cm
	0.5 mm	17.0 cm	20.4 cm	23.8 cm	27.2 cm	30.6 cm	34.0 cm
	0.6 mm	11.8 cm	14.1 cm	16.5 cm	18.9 cm	21.2 cm	23.6 cm
	0.8 mm	6.6 cm	8.0 cm	9.3 cm	10.6 cm	11.9 cm	13.3 cm
	1.0 mm	4.2 cm	5.1 cm	5.9 cm	6.8 cm	7.6 cm	8.5 cm

4.5 Precautions to Take Before Running the KOBRA® CELL

- This device has only to be used with the supplied power pack.
- Always ensure there is a mobile phase flowing through the KOBRA® CELL before switching on the current source.
- Do not switch on without mobile phase flowing through the KOBRA® CELL or the membrane will be damaged.
- It is advised to run the mobile phase through the KOBRA® CELL for 30 minutes prior to analysis.
- Always turn off the KOBRA® CELL power pack before turning off the HPLC pump.
- Always use HPLC grade solvents (with minimal benzene content) from a quality supplier and maintain fresh supplies of potassium bromide.
- Never flush 100 % solvent through the KOBRA® CELL as this can damage the membrane.

4.6 Operating the KOBRA® CELL

The KOBRA® CELL power pack is turned on by moving the switch to the 'On' position. To shut down the device move the switch to the 'Off' position.

4.7 Error LED

- The LED will light up if there is no KOBRA® CELL attached to the current source.
- The LED will also light up when the KOBRA® CELL is dry and the flow of the current source is inhibited.

Power Light

Error Light



5 HPLC Conditions

HPLC Conditions	
Derivatisation	KOBRA® CELL at 100 µA setting
Guard Cartridge	Inertsil ODS-3 5 µm, 4 mm x 10 mm or (Hichrom) equivalent
Analytical Column	Inertsil ODS-3V 5 µm, 4.6 mm x 150 mm (Hichrom) or equivalent
Mobile Phase	Water : Methanol (60 : 40 v/v) Add 119 mg of potassium bromide and 350 µl 4 M Nitric Acid to 1 litre of mobile phase
HPLC Pump	To deliver mobile phase
Flow Rate	1.0 ml/minute
Fluorescence Detector	Excitation: 362 nm Emission: 425 nm (B1 and B2) 455 nm (G1 and G2)
Column Heater	Maintain guard and analytical columns at 40 °C
Integrator / Data Control System	From preferred supplier
Injector	Autosampler / Reodyne valve
Injection Volume	100 µl
Elution Order	G2, G1, B2, B1.

6 KOBRA® CELL Maintenance

6.1 Storage

- The KOBRA® CELL should be switched off and disconnected for washing and storage (Fig. 5).
- Always store the KOBRA® CELL filled with water in order to keep the membrane wet.

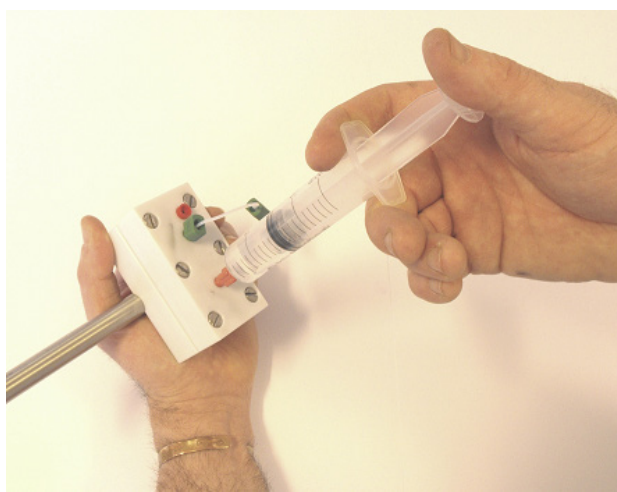


Fig 4. Storage of KOBRA® CELL

To fill the KOBRA® CELL with water:

1. Fill a syringe with 5 - 10 ml of water and attach the tip of the syringe to a pipette tip.
2. Unscrew one of the storage plugs / ferrules from one inlet / outlet position of the KOBRA® CELL, leaving the other inlet / outlet positions sealed with storage plugs or ferrules.
3. Using the syringe with the pipette tip, insert this into the empty inlet / outlet position on the KOBRA® CELL and fill the cell with water until water overflows from the syringe connection.
4. When finished, reinsert the storage plug / ferrule to contain the water within the cell and keep the membrane wet.

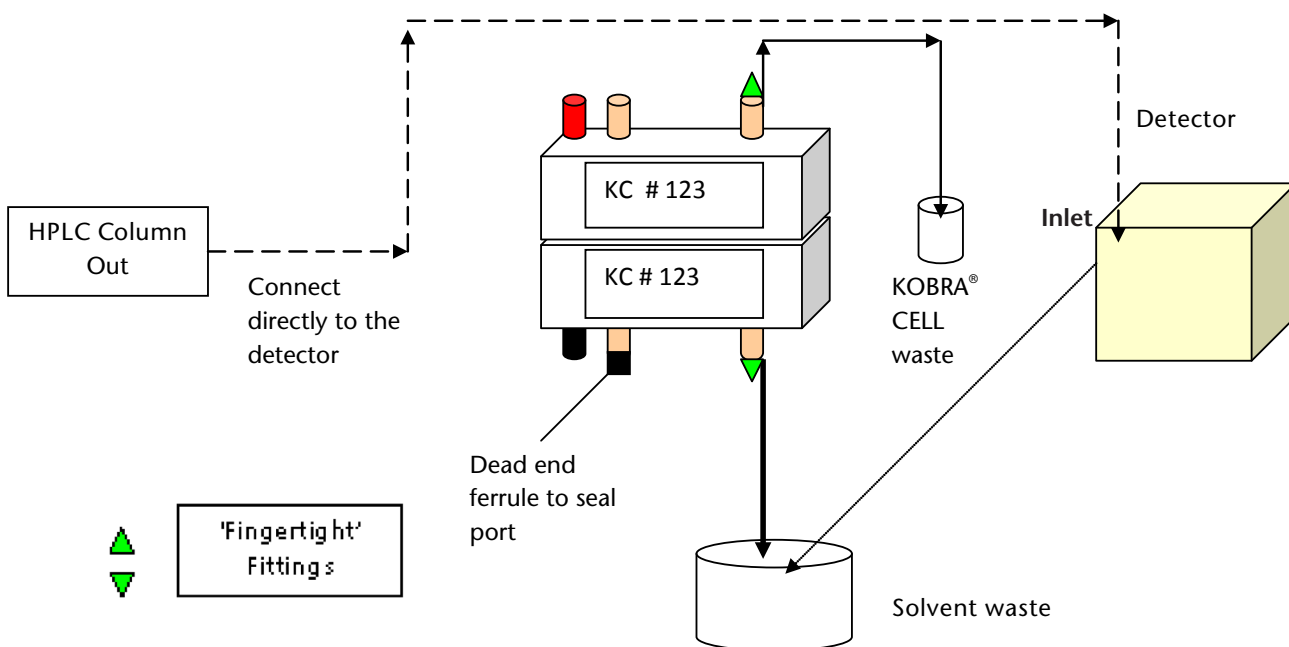


Fig 5. KOBRA® CELL Switched Off for Washing and Storage

6.2 Daily Cleaning

The HPLC system can be left running overnight at a reduced flow rate (e.g. 0.1 ml/min) and with the fluorescence detector and KOBRA® CELL power source turned off, if necessary. Alternatively, a better practice to prolong the life of the system is to clean each day, as follows:

1. Turn off KOBRA® CELL power pack and fluorescence detector.
2. Switch off HPLC pump.
3. Disconnect the KOBRA® CELL from the HPLC system and reconnect the HPLC column directly to the fluorescence detector.
4. Change the mobile phase to 5 % solvent and put the HPLC column discharge to waste.
5. Switch on the HPLC pump and flush the system for at least 30 minutes before switching off the system.
6. Manually flush the KOBRA® CELL through with distilled water several times using a syringe.
7. Store the KOBRA® CELL overnight with water inside the cell by closing off both ports on each housing using the storage plugs and capillary tubing supplied with the cell.
8. The next day, change the mobile phase from 5 % solvent back to the normal mobile phase of potassium bromide and nitric acid, and flush the HPLC system through for at least 30 minutes before re-connecting the KOBRA® CELL.
9. Reconnect the KOBRA® CELL (Fig. 6).

Caution: 100 % acetonitrile will damage the KOBRA® CELL membrane.

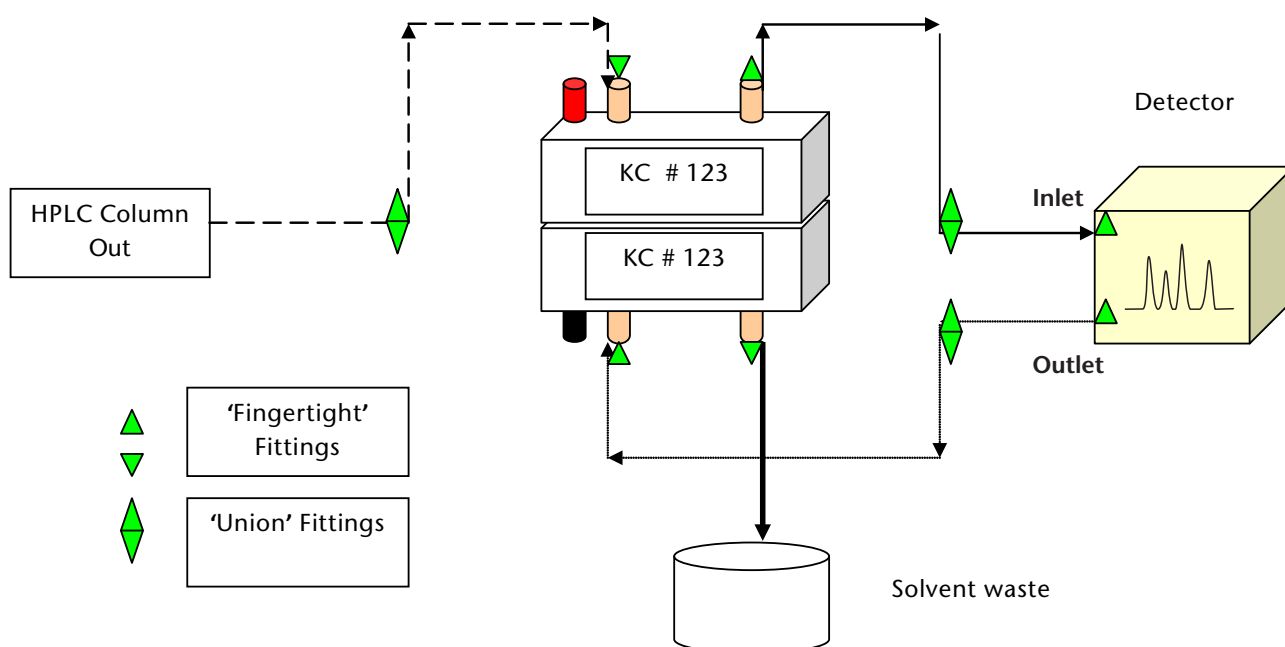


Fig 6. KOBRA® CELL Switched On for Analysis

6.3 Monitoring Performance

It is necessary to regularly monitor the performance of the KOBRA® CELL in order to detect any deterioration in the membrane contained within. The performance should be checked at the time of installation and then weekly by comparing the peak areas of a known aflatoxin standard. The same should also be done in order to monitor the deterioration of the lamp in the detector.

Deterioration of the membrane will occur over a period of time depending upon the frequency of use and the type of samples analysed (Fig. 7). When the performance becomes unacceptable the membrane should be replaced (Product Code: K02).

Normally, it is found that even under extreme workloads the membrane does not need to be replaced for at least 6 months or 1,000 injections.



Fig 7. Deterioration of KOBRA® CELL Membrane

6.4 Changing the Membrane

Change the membrane when a decrease in performance is observed which cannot be attributed to the detector lamp. Additional membranes are available separately (K02).

1. Using a screw driver take off the white plastic housing of the KOBRA® CELL by removing the 6 locating screws.
 2. Carefully separate the top white plastic housing, and remove in turn the internal layers of the cell using tweezers. Make a note of the position and orientation of each layer as it is removed.
 3. Continue removing the internal layers until the membrane is exposed. Before removal inspect the membrane for damage, discolouration and orientation. During the life of the KOBRA® CELL the bromination reaction causes the centre of the membrane to lose it's colour. No holes should be visible in the membrane other than the pre-cut holes for the locating pins at either end of the membrane. Note the orientation of the membrane.
 4. Remove the membrane from the KOBRA® CELL ready for replacement using tweezers.
 5. Carefully remove the spare membrane, holding it by one end with the tweezers rather than the middle.
 6. Position the replacement membrane over the different sized locating pegs, according to the size of pre-cut holes in the membrane.
- Caution:** Do not allow the membrane to dry out. Add distilled water if necessary.
7. Carefully replace the layers to be positioned on top of the membrane, remembering their correct order and orientation.
 8. Secure the layers using the 6 locating screws.

7 Troubleshooting

7.1 No Peaks Observed

- If no peaks are observed on the HPLC chromatogram it is probably not the KOBRA® CELL that is at fault (Fig. 8). There may be a problem with a component of the HPLC system or with the standards used.
- Check that the bulb in the lamp does not need replaced. It should be replaced every 1,000 - 2,000 hours depending on the detector.
- Check there is not a problem with the detector.
- Check the standard has been prepared correctly, or the standard has not deteriorated. Standards should be prepared fresh every day.
- Check the PEEK™ tubing leading to and from the KOBRA® CELL is not blocked or distorted. Replace if required.
- Check the internal components of the KOBRA® CELL have been inserted in the correct order, particularly if no peaks are observed following the replacement of one of the components (Fig. 1).

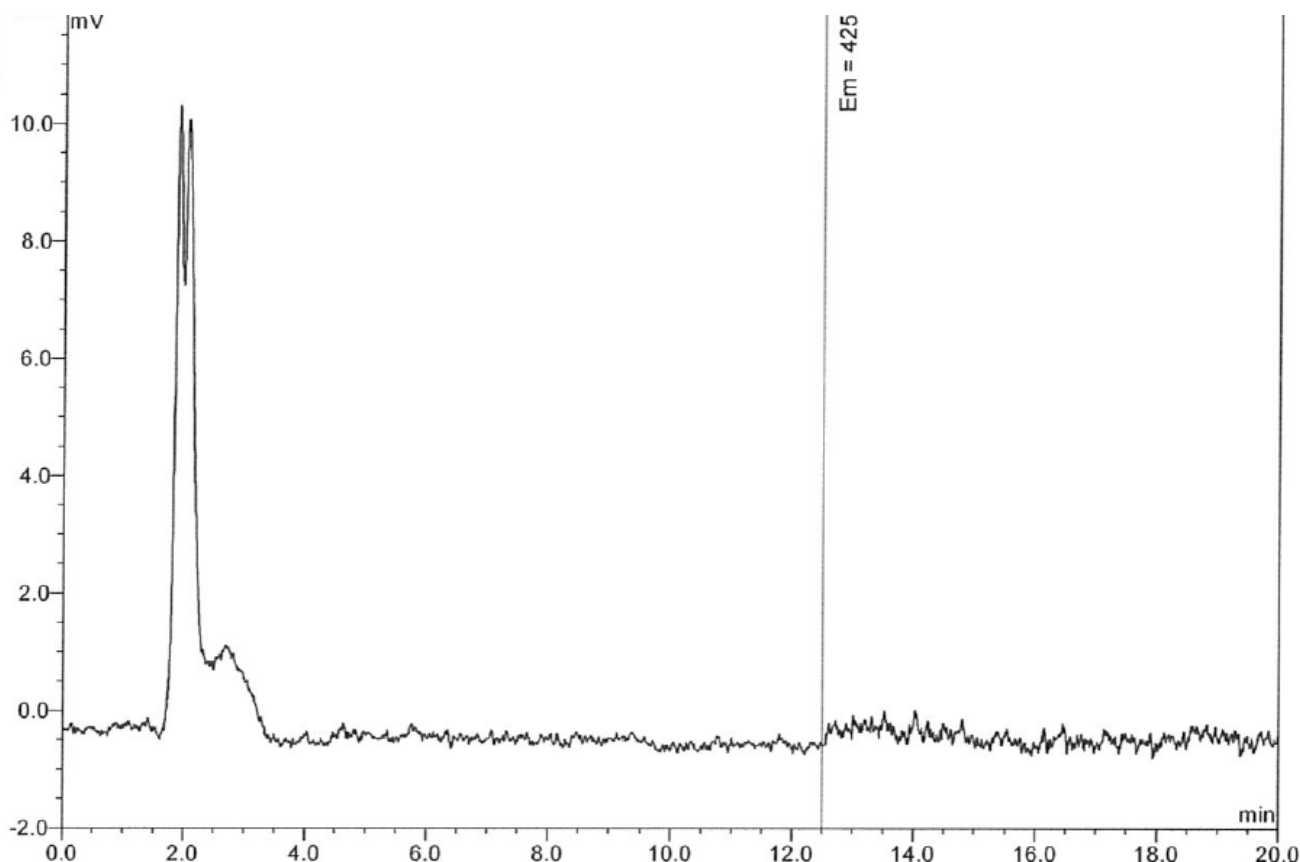


Fig 8. No Peaks Obtained in Chromatogram

7.2 Reduced / No Aflatoxin B1 and G1 Peaks Observed

- Check the KOBRA® CELL is connected correctly to the HPLC (Fig 6).
- Check the connections to the KOBRA® CELL from the power pack.
- Check the KOBRA® CELL membrane. The membrane is the site of the actual derivatisation reaction (ion-exchange occurs on the membrane). The membrane starts a deep amber colour and, through time and use, the centre of the membrane becomes discoloured, eventually turning a pale colour. The pale colour indicates that the membrane will need replaced.
- Check there is power in the power pack. Check that the power pack has been turned on.
- Check the internal components of the KOBRA® CELL are in the correct order (Fig. 1), or that there is not a blockage in the system.
- No potassium bromide (or not enough potassium bromide) was added to the mobile phase. Prepare again.
- No nitric acid was added to the mobile phase. Prepare again.

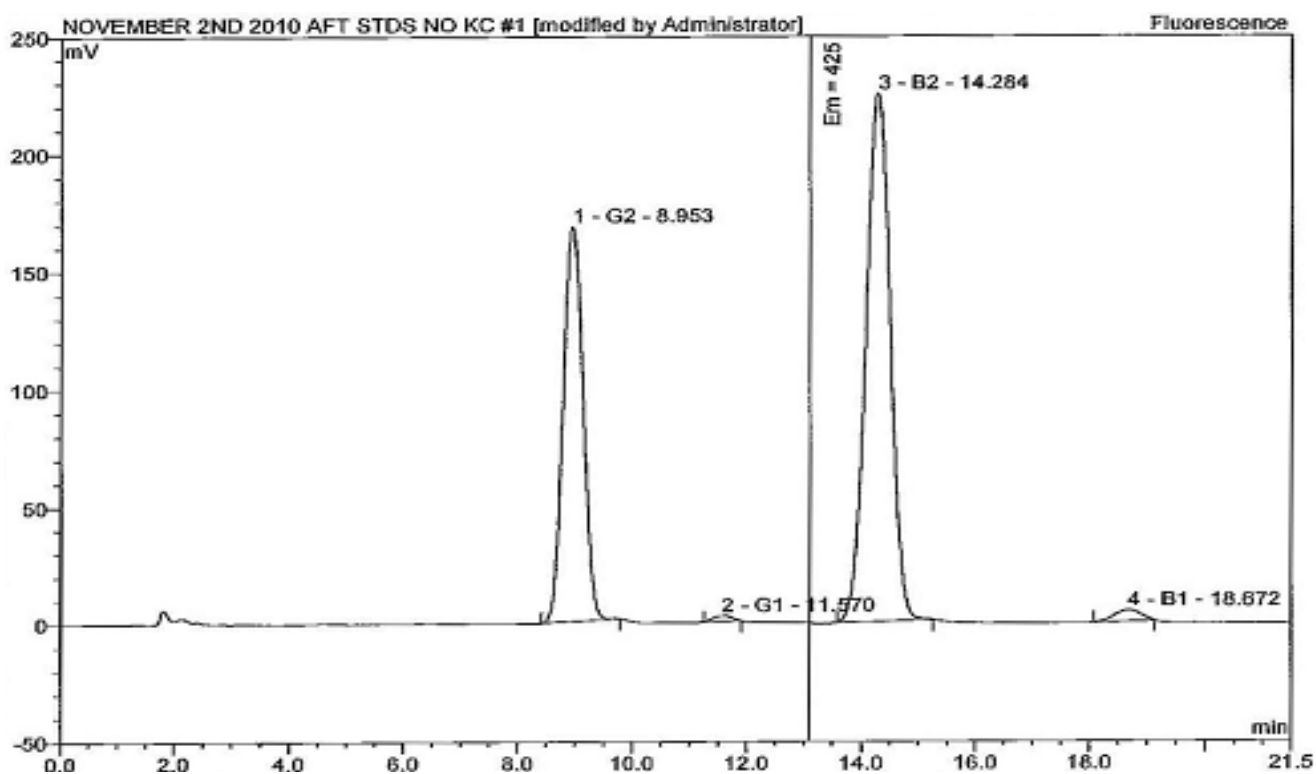


Fig 9. Reduced / No Aflatoxin B1 and G1 Peaks Observed

7.3 Abnormal Peaks Observed

- If abnormal peaks are observed on the HPLC chromatogram it is not normally the KOBRA® CELL that is at fault (Fig. 10).
- Check the correct mobile phase has been used.
- Check the correct HPLC column has been used.
- Check the internal components of the KOBRA® CELL have been inserted in the correct order (Fig. 1).
- Analytical column could need changed and may result in poor peak resolution or peak tailing.

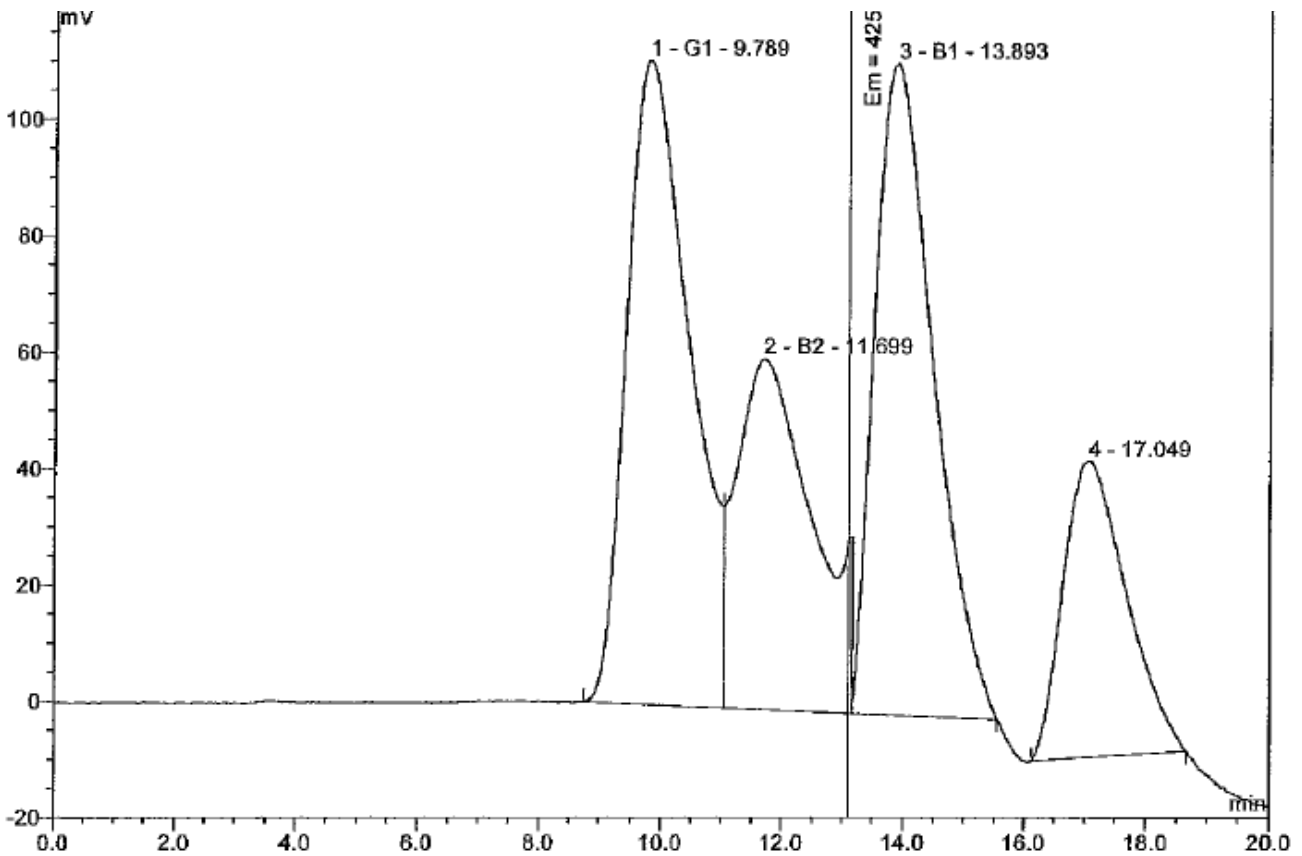


Fig 10. Abnormal Peaks Obtained in Chromatogram

7.4 Error Light Observed on Power Pack

- The LED will light up if there is no KOBRA® CELL attached to the current source.
- The LED will also light up when the KOBRA® CELL is dry as the flow of the current source is inhibited.
- Check the red and black terminal leads are connected to the KOBRA® CELL correctly.
- Check the KOBRA® CELL has mobile phase running through it. It is important to keep liquid in the KOBRA® CELL at all times in order for it to function. It is advised to run mobile phase through the KOBRA® CELL for at least 30 minutes prior to injection of the standards or samples.
- Check the power pack. If a connection or component of the power pack is faulty the power pack should be replaced.

7.5 Leakage from KOBRA® CELL

- The locating screws on the KOBRA® CELL have not been tightened correctly. Please refer to 3.2, Precautions When Assembling.

- Check the locating screws on the KOBRA® CELL have not been over tightened causing damage to the threads. Once the threads are damaged a new white plastic housing is required.
- Check the PEEK™ tubing is not blocked or distorted. Replace if required.

7.6 KOBRA® CELL Causing High Pressure on the HPLC

- Check the PEEK™ tubing is not blocked or distorted. Replace if required.
- Check the inlet and outlet connections that the PEEK™ tubing is inserted into. If they are blocked they can be unscrewed and sonicated.
- Check the internal components in the KOBRA® CELL have been inserted in the correct orientation. Reassemble if required (Fig. 1 and 2).

7.7 Storage Plug Broken Off or Lodged in KOBRA® CELL

- Return to R-Biopharm Rhône for repair.

8 General Information

8.1 Quality

RBR products are developed, manufactured, tested and dispatched under an ISO 9001 and ISO 13485 registered Quality Management System, guaranteeing a consistent product, which always meets our performance specifications. Our products have been used in many collaborative studies to develop standard European and International Methods and are widely used by key institutions, food companies and government laboratories. Customer references for RBR products are available on request.

8.2 Technical Support

RBR understand that from time to time users of our products may need assistance or advice. Therefore, we are pleased to offer the following services to our customers:

- Analysis of problem samples.
- Application notes for difficult samples.
- References from the RBR library.
- Installation and support of the KOBRA® CELL.
- Advice on detection parameters.
- Advice on preparation and handling of standards.
- Updates on legislation, sampling and other news by e-mail.
- Provision of spiked samples.

Please contact your local R-Biopharm distributor for further information.

8.3 Warranty

R-Biopharm Rhône Ltd makes no warranty of any kind, express or implied, except that all products made by R-Biopharm Rhône Ltd are made with materials of suitable quality. If any materials are defective, R-Biopharm Rhône Ltd will provide a replacement product. The user assumes all risk and liability resulting from the use of R-Biopharm Rhône Ltd products and procedures. R-Biopharm Rhône Ltd shall not be liable for any damages, including special or consequential damages, loss or expense arising directly or indirectly from the use of R-Biopharm Rhône Ltd products or procedures.

