

STREPTAVIDIN

SPECIFICATIONS

Product Code: SA10
Specific Activity: ≥ 14.0 U/mg
Purity: Predominantly a single band by SDS-PAGE.
Shipped with ice pack for next day delivery.
Storage: -20°C
Formulation: Lyophilized (contains ~0.9 mg protein/mg lyophilizate; balance is sodium chloride).
Stability: Lyophilized streptavidin is stable for at least one year when stored desiccated at -20°C .

Streptavidin (CAS No. 9013-20-1) is a biotin-binding protein found in the culture broth of the bacterium *Streptomyces avidinii*. Like its namesake avidin, streptavidin binds 4 moles of biotin per mole of protein with a high affinity virtually unmatched in nature ($K_d \sim 10^{-15}$). Streptavidin lacks carbohydrate side chains present on avidin and has an isoelectric point nearer to neutrality where most useful biological interactions occur (pI of 5-6 *vs* 10 for avidin). As a result, streptavidin frequently exhibits lower levels of non-specific binding than does avidin when the proteins are used in applications relying upon the formation of avidin/biotin complexes (Chalet and Wolf, 1964).

In ELISA-based diagnostic systems, antibodies directed against a particular antigen may be covalently attached to

reporter enzymes. Antigens are then quantitated by enzymatic assay after binding to these conjugated molecules.

Unfortunately, the precise conditions for accomplishing such covalent attachments must be determined individually for each antibody/reporter combination, and often result in significant loss of either the enzymatic activity of the reporter enzyme or the binding functions of the antibodies.

Streptavidin finds utility in these systems because antibody molecules are easily modified by the covalent attachment of derivatives of biotin with little or no loss in the ability of the antibody molecules to bind their antigens. These biotinylated antibodies may be detected by their interaction with conjugates of streptavidin and the reporter enzymes (Guesdon *et al.*, 1979). The same preparation of conjugated streptavidin-reporter enzyme may be used with any number of different biotinylated antibodies making this system a highly flexible one.

The reporter molecule may be bound to streptavidin covalently, or biotinylated and attached to streptavidin via the streptavidin-biotin interaction. Since streptavidin is multivalent (binding 4 molecules of biotin per tetrameric protein molecule) it may be used in combination with biotinylated antibody and biotinylated reporter enzymes to obtain amplified signals. Such amplification in ELISA's is otherwise difficult to obtain and requires the introduction of additional antibody components. ELISA systems employing streptavidin can readily detect sub-nanogram amounts of antigens.

CHARACTERISTICS

Molecular Weight: 52,000 daltons, composed of 4 essentially identical polypeptide chains.

Composition: Streptavidin is a homotetramer with 24-32 lysines per tetramer. The protein contains no cysteine residues, carbohydrate side chains or associated cofactors. Different preparations of streptavidin show considerable heterogeneity at both the amino- and carboxy-termini of each subunit polypeptide due to proteolysis during biosynthesis and secretion. Monomeric subunits of streptavidin are synthesized as 183 amino acid prepeptides. During secretion by *Streptomyces* sp., a 24 amino acid leader sequence is cleaved from these polypeptides resulting in newly secreted monomers of 159 amino acids (Agaraña *et al*, 1986). Upon longer incubations in culture these monomers are progressively cleaved to "core" subunits containing 125-127 amino acids (Bayer *et al*, 1989; Pähler *et al*, 1987). Preparations of streptavidin are relatively stable over a wide pH range and extremely heat stable, requiring up to 20 minutes at 100°C in 0.2% SDS to dissociate the subunits (Bayer *et al*, 1986). Strong chaotropic agents such as 6 *M* urea have been reported to dissociate the streptavidin tetramer into dimers (Sano and Cantor, 1990). These dimers appear to be stable in urea without appearance of monomers.

Binding Properties: Unproteolyzed and proteolyzed preparations of streptavidin appear to bind biotin with equal affinity. The most highly proteolyzed tetramers may bind over 16 micrograms of d-biotin per milligram of protein. Bayer (1989) reports that biotinylated enzymes bind most effectively to truncated streptavidin in ELISA-type assays.

The dissociation constant for biotin is approximately 10^{-15} *M* (Chaiet *et al*, 1964). The formation of the streptavidin-biotin complex is stable over wide pH and temperature ranges. The complex is generally disrupted only by conditions which lead to irreversible denaturation of the protein. Analogs of biotin such as 2-imino-biotin bind reversibly to the protein with complex formation at high pH (>9.5) and dissociation at low pH (<4) (Bayer *et al*, 1986; Hofmann *et al*, 1980).

Extinction Coefficient: $E_{280}^{1\%} = 32$ (Suter *et al*, 1988)

Homogeneity: The streptavidin gene product is subject to cleavage by endogenous proteases during fermentation. In some instances this can result in a heterogeneous product. The crude material may be treated with exogenous proteases during purification to reduce this heterogeneity. ProZyme Streptavidin is a pure, homogeneous preparation which shows predominantly a single band by SDS PAGE.

Applications: Enzyme conjugates and complexes; Southern blots and other methodologies related to DNA and RNA analysis; Western blots; and purification of proteins or other antigens with biotinylated antibodies or lectins by use of immobilized streptavidin.

Origin: USA

ASSAY

The biotin-binding activity of streptavidin is determined using a modification of the dye-binding assay of Green (1970). One unit will bind one microgram of d-biotin at pH 7.0.

Reagents

10 mM 2-(4'-hydroxyazobenzene)benzoic acid (Sigma) dissolved in 10 mM sodium hydroxide (HABA)

0.2 M sodium phosphate, pH 7.0

2 mM d-biotin in 0.1 M sodium phosphate, pH 7.0

Streptavidin sample - Dissolve at 5-10 mg/ml in de-ionized water.

Procedure

Adjust spectrophotometer to read at 500 nm.

To two tubes labeled A and B add as follows:

	A	B
Streptavidin sample	0.05 ml	0.05 ml
Phosphate Buffer	0.5 ml	0.5 ml
HABA stock	0.1 ml	0.1 ml
Biotin stock	--	0.25 ml
H ₂ O	0.35 ml	0.1 ml
Total Volume	1.0 ml	1.0 ml

After mixing, zero the spectrophotometer with water and read the absorbance values in tubes A and B.

Calculation

$$\text{Units/mg} = \frac{(10^6 \mu\text{g/g})(A-B)MV}{\epsilon C v} = \frac{141(A-B)}{C}$$

where:

M = formula weight of d-biotin (244 g/mole)

V = volume of assay in liters (.001 liters)


v = volume of streptavidin sample in milliliters (.05 ml as written)

C = concentration of streptavidin in sample (mg/ml)

ϵ = net molar extinction coefficient of HABA-streptavidin complex at 500 nm (34,500 M^{-1})

The measured specific activity of a given lot depends on the assay method used, as well as the method used to quantitate the Streptavidin (gravimetric *vs.* spectrophotometric). For comparison purposes, ProZyme also reports the results of a biotin titration assay.

Biotin-binding capacity:

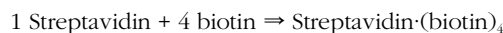

The formation of the complex is measured at 233 nm.

Reagents

Ammonium Carbonate Buffer:

200 mM $(\text{NH}_4)_2\text{CO}_3$, pH 8.9

1 mg/ml solution of d-biotin (dissolved in Ammonium Carbonate Buffer)

10 mg/ml Streptavidin sample (dissolved in de-ionized water)

Procedure

Adjust spectrophotometer to read at 233 nm.

Dilute the dissolved streptavidin sample to 0.1 mg/ml in Ammonium Carbonate Buffer. Prepare at least 1 ml.

Dilute the 1 mg/ml biotin solution to ~0.1 mg/ml in Ammonium Carbonate Buffer. Prepare at least 1 ml.

Add 1 ml of 0.1 mg/ml Streptavidin sample to cuvette. Take initial reading of unbound streptavidin.

Add diluted biotin (0.1 mg/ml) in increments of 3 μl . Mix carefully after each addition and take absorbance readings.

Continue adding biotin until absorbance readings plateau or decline repeatedly.

Graph titration curve as absorbance (ordinate) *vs.* volume of biotin added (abscissa). Note the volume (μl) of biotin added at the inflection point of the curve.

The amount of biotin at the inflection point is divided by the concentration of the streptavidin sample to yield the specific activity.

Calculation

$$\text{Units/mg} = \frac{V_1 \times C_1}{V_2 \times C_2}$$

where:

V_1 = vol of d-biotin at inflection point

C_1 = concentration of d-biotin solution

V_2 = volume of streptavidin sample

C_2 = concentration of streptavidin sample

SUGGESTIONS FOR USE

Measurement

Because the lyophilizate is ~10% NaCl by weight, we suggest that the customer dissolve it and measure it using its extinction coefficient [$E_{280}(1\%) = 32$] rather than weighing it.

Aliquots may be stored at -20°C for long periods. If kept refrigerated, an antimicrobial agent such as 0.05% sodium azide should be added to retard growth and subsequent proteolysis from exogenous enzymes.

Solubility

ProZyme Streptavidin is readily soluble in water or salt-containing buffers, up to 50 mg/ml or more.

Bayer (1989) reports that streptavidin may form aggregates under certain conditions. There is a tendency for lyophilized streptavidin to aggregate when it is re-dissolved in water or other low ionic strength buffers at neutral or acidic pH.

ProZyme Streptavidin has been lyophilized from a dilute sodium chloride solution at mildly alkaline pH in order to minimize aggregate formation.

However, in rare cases Streptavidin may contain a small amount of insoluble material when dissolved in de-ionized water or low ionic strength buffers, either at the time it is dissolved or after freezing and re-thawing. This effect is generally not seen in the presence of salt-containing buffers such as PBS.

If undissolved material is observed, it can be removed by centrifugation and does not constitute a significant fraction of the total protein. The activity of the material recovered after reconstitution under these conditions is undiminished.

REFERENCES

- Agaraña, C. E., Kuntz, I. D., Birken, S., Axel, R. and C. R. Cantor. Molecular cloning and nucleotide sequence of the streptavidin gene. **Nucleic Acids Res** **14**: 871-882 (1986).
- Bayer, E. A., Ben-Hur, H., Hiller, Y. and M. Wilchek. Postsecretory modifications of streptavidin. **Biochem J** **259**: 369-376 (1989).
- Bayer, E. A., Ben-Hur, H., Gitlin, G. and M. Wilchek. An improved method for the single-step purification of streptavidin. **J Biochem Biophys Meth** **13**: 103-112 (1986).
- Brower, M. S., Brakel, C. L. and K. Garry. Immunodetection of streptavidin-acid phosphatase complex on western blots. **Analyt Biochem** **147**: 382-386 (1985).
- Cartun, R. W. and C. A. Pedersen. An immunocytochemical technique offering increased sensitivity and lowered cost with a streptavidin-horseradish peroxidase conjugate. **J Histotechnology** **12**: 273-277 (1989).

- Chaiet, L. and E. J. Wolf. The properties of streptavidin, a biotin-binding protein produced by *Streptomyces*. **Arch Biochem Biophys** **108**: 1-5 (1964).
- Gitlin, G., Bayer, E. A. and M. Wilchek. Studies on the biotin-binding site of streptavidin. **Biochem J** **256**: 279-282 (1988).
- Green, N. M. Spectrophotometric determination of avidin and biotin, in **Meth Enzymol** **18** (McCormick, D. B. and L. D. Wright, eds.) Academic Press, New York, p. 418 (1970).
- Green, N. M. Avidin. **Adv Protein Chem** **29**: 85-133 (1975).
- Guesdon, J.-L., Ternynck, T. and S. Avrameas. The use of avidin-biotin interaction in immunoenzymatic techniques. **J Histochem Cytochem** **27**: 1131-1139 (1979).
- Hofmann, K., Wood, S. W., Brinton, C. C., Montibeller, J. A. and F. M. Finn. Iminobiotin affinity columns and their application to retrieval of streptavidin. **Proc Natl Acad Sci USA** **77**: 4666-4668 (1980).
- Leary, J. J., Brigati, D. J. and D. C. Ward. Rapid and sensitive colorimetric method for visualizing biotin-labeled DNA probes hybridized to DNA or RNA immobilized on nitrocellulose: Bio-blots. **Proc Natl Acad Sci USA** **80**: 4045-4049 (1983).
- Moxley, G. Heightened sensitivity of quantitative ELISA for IgM rheumatoid factor with the use of the biotin-streptavidin system. **Am J Clin Path** **92**: 630-636 (1989).
- Nikolau, B. J., Wurtele, E. S. and P. K. Stumpf. Use of streptavidin to detect biotin-containing proteins in plants. **Analyt Biochem** **149**: 448-453 (1985).
- Pähler, A., Hendrickson, W. A., Gawinowicz Kolks, M. A., Argaraña, C. E. and C. R. Cantor. Characterization and crystallization of core streptavidin. **J Bio Chem** **262**: 13933-13937 (1987).
- Sano, T. and C. R. Cantor. Cooperative biotin binding by streptavidin. **J Biol Chem** **265**: 3369-3373 (1990).
- Sharma, H. and C. Tihon. Determination of avidin and streptavidin by a modified Bradford method. **Analyt Biochem** **170**: 135-139 (1988).
- Suter, M., Cazin, J., Jr., Butler J. E. and D. M. Mock. Isolation and characterization of highly purified streptavidin obtained in a two-step purification procedure from *Streptomyces avidinii* grown in a synthetic medium. **J Immunol Meth** **113**: 83-91 (1988).
- Weber, P. C., Ohlendorf, D. H., Wendoloski, J. J. and F. R. Salem. Structural origins of high-affinity biotin binding to streptavidin. **Science** **243**: 85-88 (1989).



3832 Bay Center Place
Hayward, CA 94545-3619

TOLL FREE (800) 457-9444
PHONE (510) 638-6900
FAX (510) 638-6919

E-MAIL info@prozyme.com
WEB www.prozyme.com