

Cat.-No.	Size
AB-102S	500 mg
AB-102L	1 g
AB-102XL	5 g
AB-102XXL	10 g

**Description:**

Nourseothricin antibiotic


**Storage conditions:**

Store at 4 °C.

Stable for 1 year.

**Selection:**

For selection of recombinant *Leishmania* strains NTC is dissolved in aqua dest. (stock solution 10 mg/ml), filtersterilized and added to growth medium to a final concentration of 100 µg/ml. The stock solution is stable for 4 weeks at 4 °C.

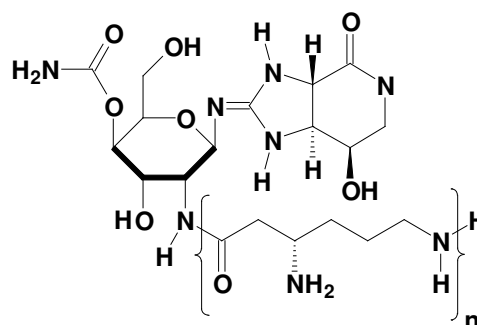

F n = 1

D n = 3

Find more information in the User's Guide for *Leishmania tarentolae* Gene Expression Starter Kit (download pdf from the LEXSY section of our website).

**FOR RESEARCH USE ONLY.**

**NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC  
OR THERAPEUTIC USES.**